


Analysis of Chlorophenols

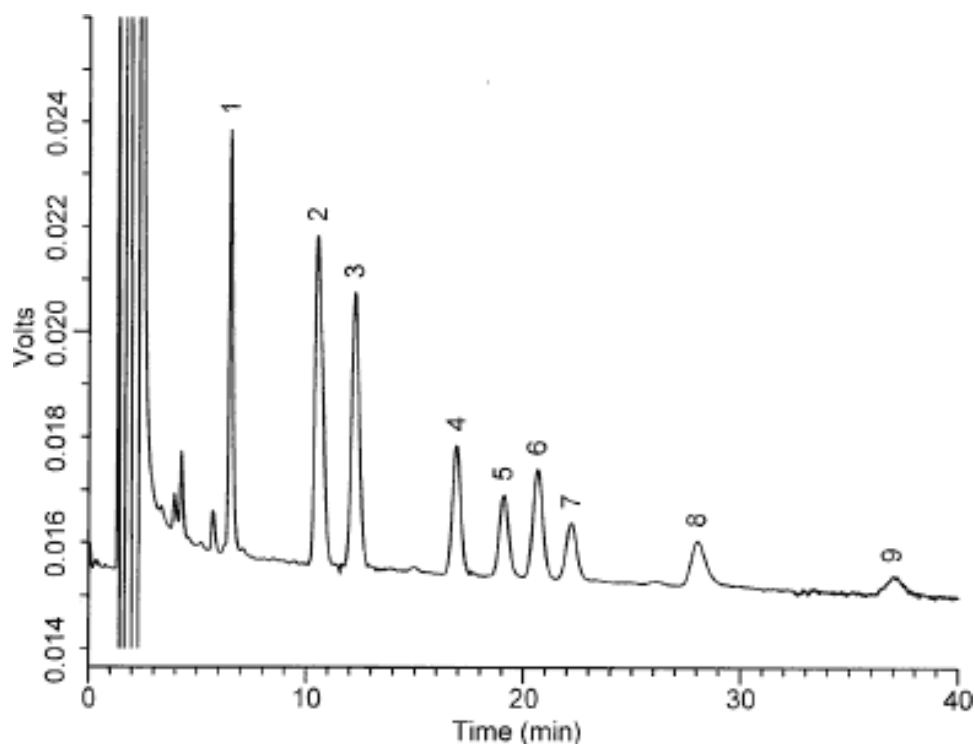

Column	: Inertsil ODS-3 150 × 1.5 mmI.D.	Solutes	: 1. 2-Chlorophenol(2-CIP)(5ng/mL) 2. 2,3-Dichlorophenol(2,3-DiCIP)(5ng/mL) 3. 2,4-Dichlorophenol(2,4-DiCIP)(5ng/mL) 4. 2,3,6-Trichlorophenol(2,3,6-TriCIP)(5ng/mL) 5. 2,3,4-Trichlorophenol(2,3,4-TriCIP)(5ng/mL) 6. 2,4,6-Trichlorophenol(2,4,6-TriCIP)(5ng/mL) 7. 2,4,5-Trichlorophenol(2,4,5-TriCIP)(5ng/mL) 8. 2,3,4,5-Tetrachlorophenol(2,3,4,5-TetraCIP)(5ng/mL) 9. Pentachlorophenol(PCP)(5ng/mL)
Column Cat. No.	: 5020-13241		
Eluent	: A) CH ₃ CN B) 50mM K ₂ HPO ₄ (pH5.0; H ₃ PO ₄) A/B = 40/60		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: ECD 750mV vs. Ag/AgCl 1000nA/V		
Injection Vol.	: 20µL		
Sample	:		