


Analysis of Hydrocarbons C2-C4 and Alcohol C1-C2

Specification

One GC Channel

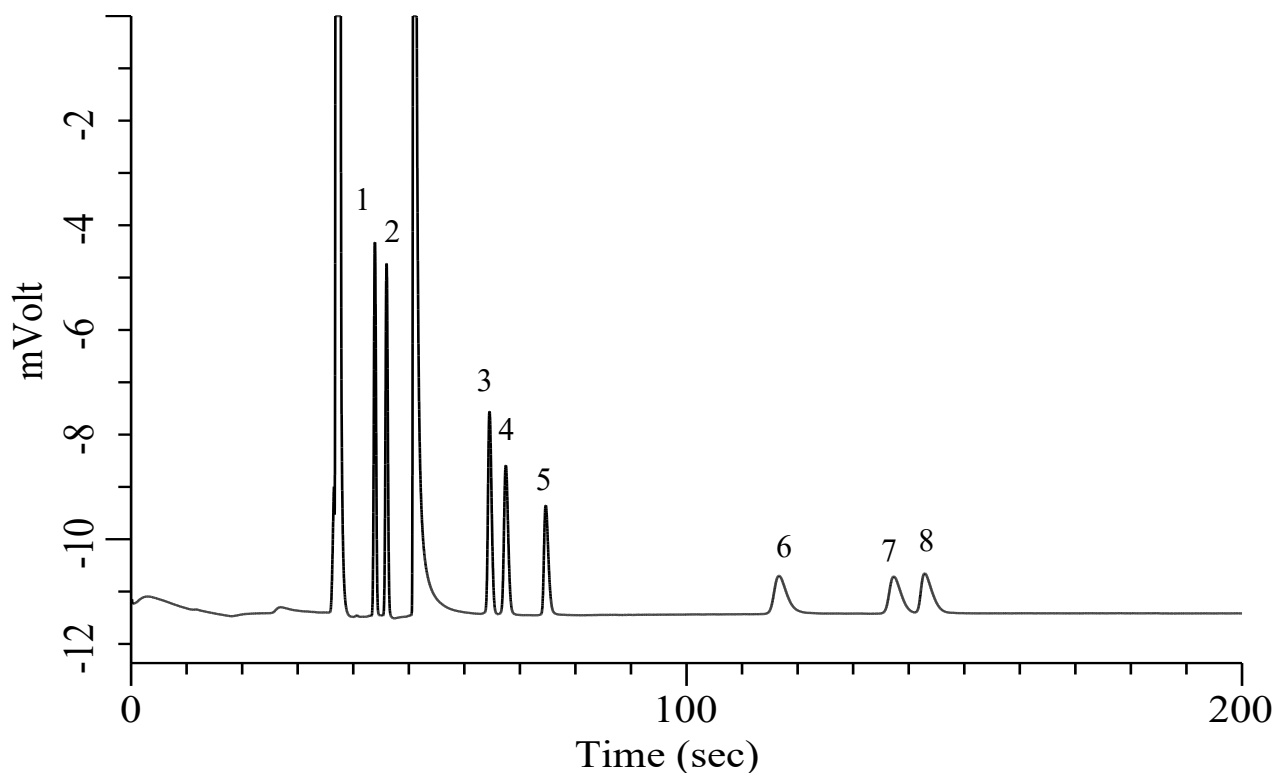

Conditions

System : Micro GC
Column : PoraPLOT Q 10 m + BF 1 m
Column Temp. : 120 °C
Carrier Gas : Ar 190 kPa
Inject Time : 40 msec

Analyte:

1. Ethylene (C_2H_4)	0.5 % (v/v)
2. Ethane (C_2H_6)	0.5 % (v/v)
3. Propylene (C_3H_6)	0.5 % (v/v)
4. Propane (C_3H_8)	0.5 % (v/v)
5. Methanol	0.5 % (v/v)
6. Isobutane ($i-C_4H_{10}$)	0.5 % (v/v)
7. <i>n</i> -Butane ($n-C_4H_{10}$)	0.5 % (v/v)
8. Ethanol	0.5 % (v/v)

Balance Gas N_2