


Vitamin E and Vitamin K

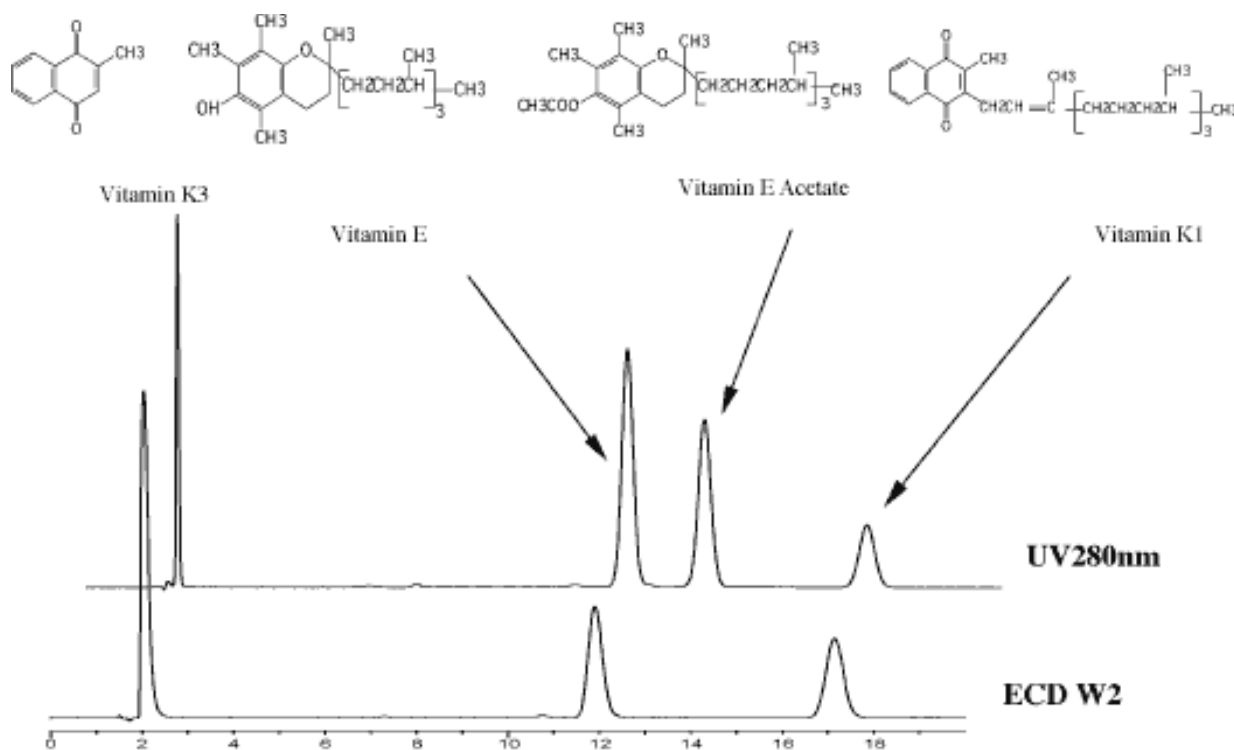

Column	: Inertsil ODS-3 5 μ m 150 $\times$ 4.6 mm I.D.	Solutes	: 1. Vitamin K3 2. Vitamin E 3. Vitamin E acetate 4. Vitamin K1
Column Cat. No. :			
Eluent	: CH ₃ CN/CH ₃ OH = 75/25, 7 g/L NaClO ₄		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 280 nm ECD W1 = -500 mV, W2 = 400 mV vs. Ag/AgCl 1000 nA/V		
Injection Vol.	: 5 μ L		
Sample	: Fat-soluble Vitamins		