


Recycle preparative purification of Steroids

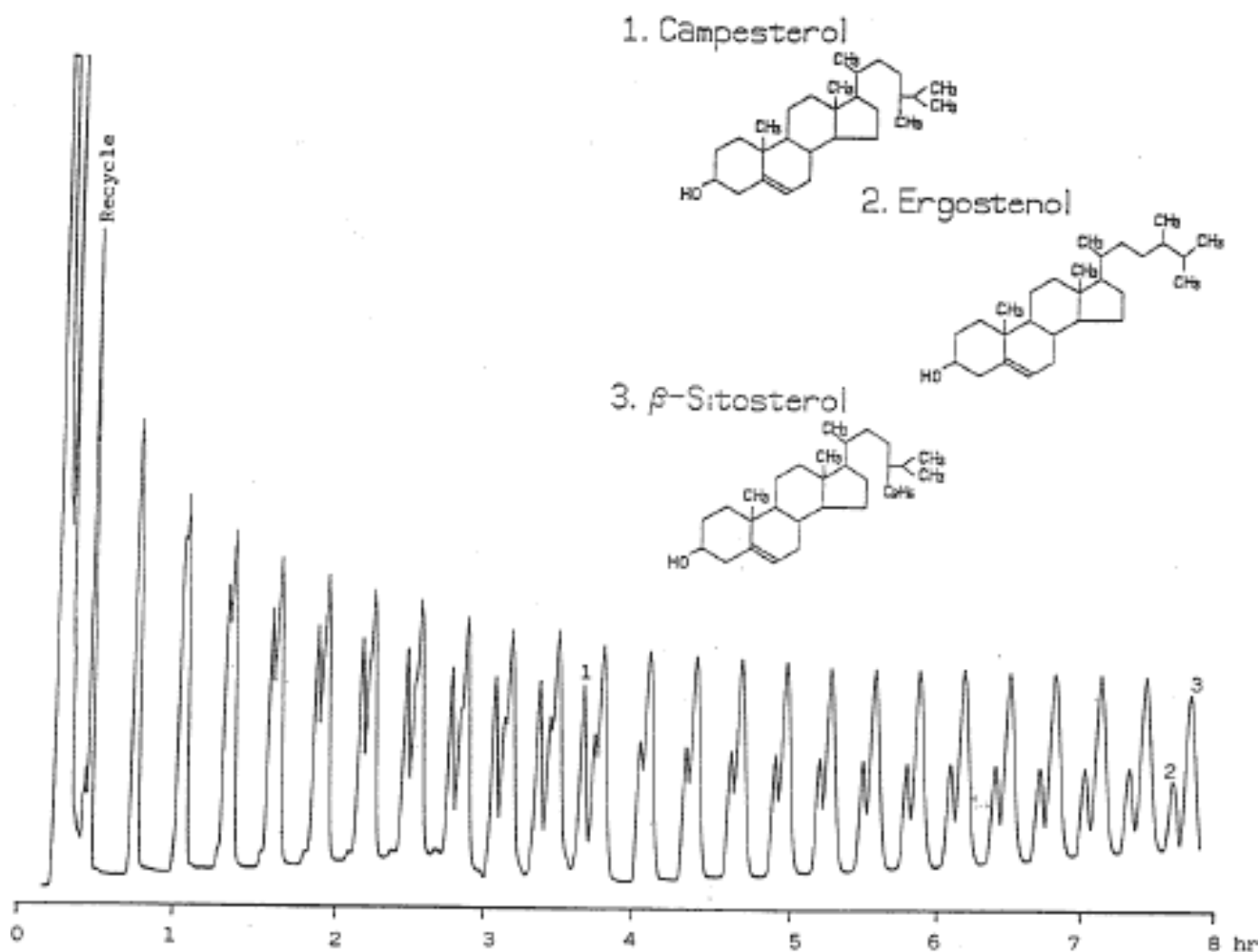

Column : Inertsil PREP-ODS

10 μ m

250 $\times$ 20mmI.D.

Solutes

:1. Campesterol

2. Ergosterol

3. β -Sitosterol

Column Cat. No. : 5020-31123

Eluent : C₂H₅OH/H₂O=93/7

Flow rate : 12 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 2mL

Sample : Crude β -sitosterol(20mg/mL)