


Thiamine

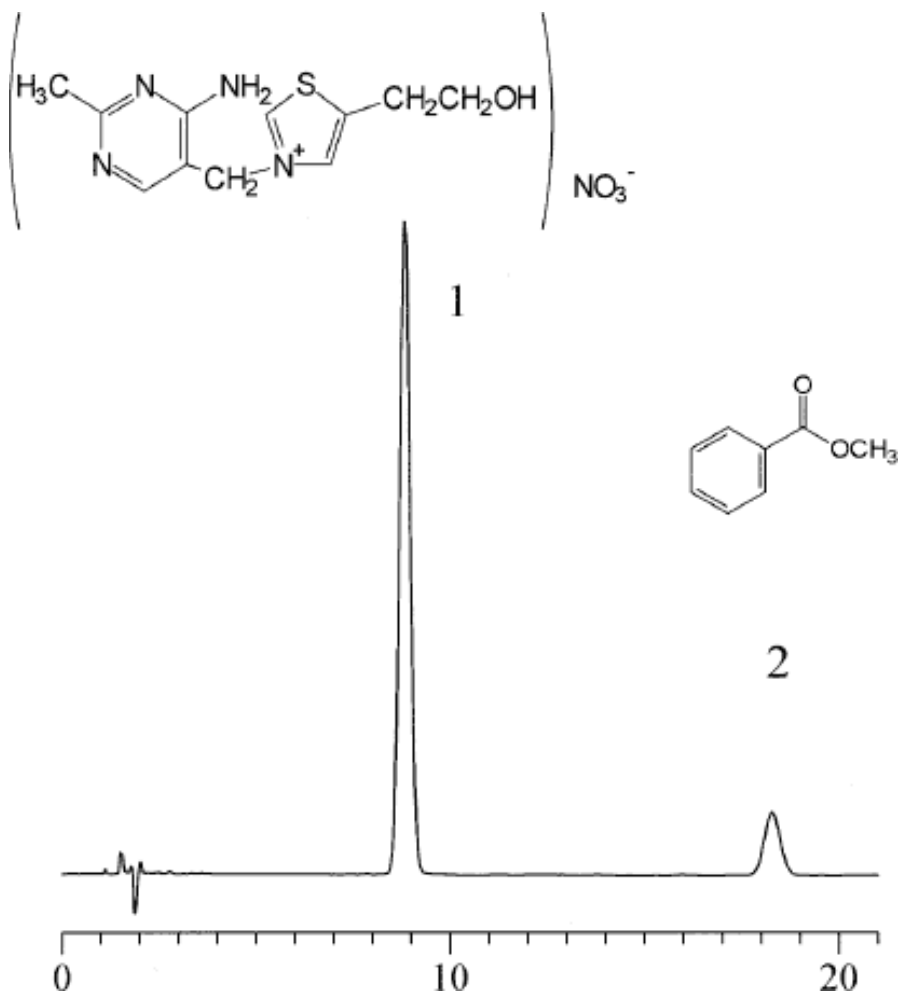

Column : Inertsil ODS-3
5 μ m
150 $\times$ 4.6mm I.D.

Column Cat. No. :

Eluent : MeOH / MeCN / 5mM Sodium 1-octanesulfonate in 1%
Acetic acid sol.
= 600 / 240 / 160

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 254 nm

Injection Vol. :

Sample : 1. Thiamine nitrate (2.0 μ g)
2. Methyl benzoate (0.4 nL)

Solutes

: 1. Thiamine nitrate (2.0 μ g)
2. Methyl benzoate (0.4 nL)