


Analysis of Oligonucleotides

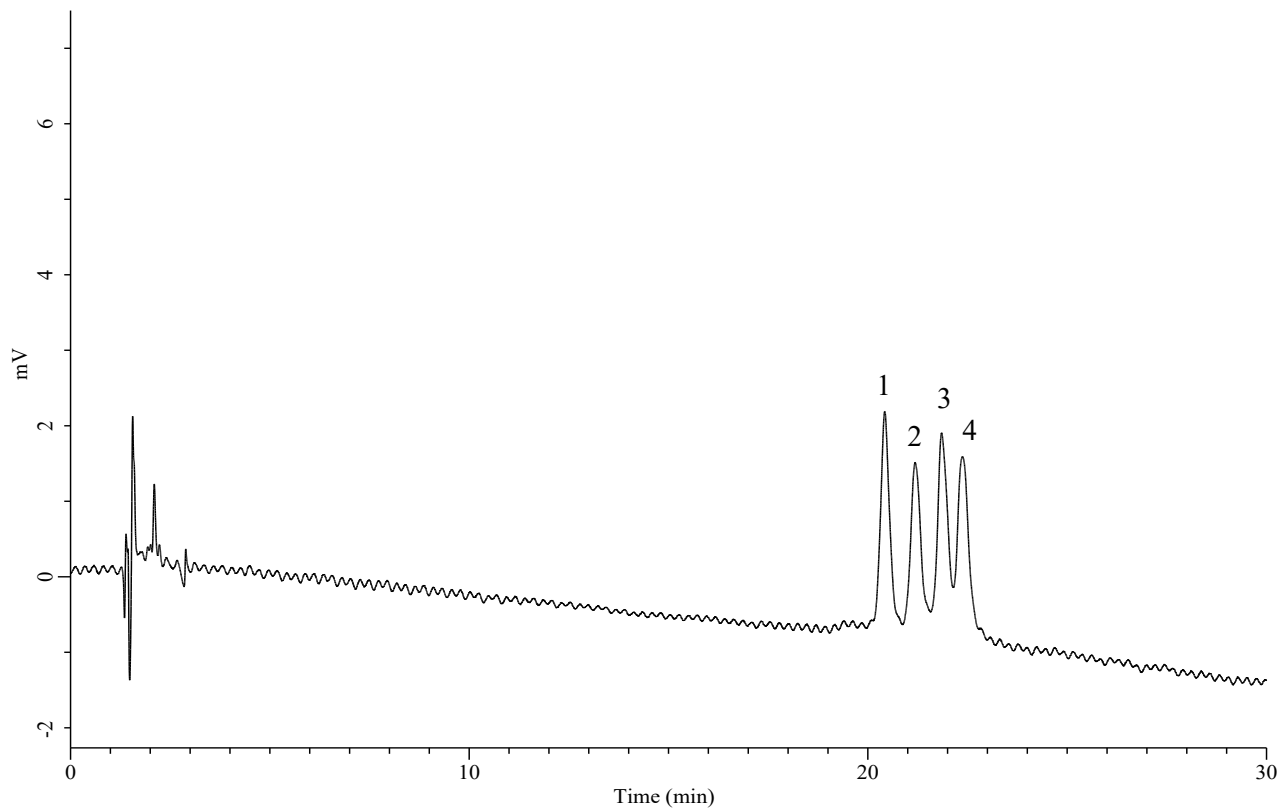

Conditions

Column : Inertsil Hybrid-Z HILIC
 (3.5 μ m, 150 $\times$ 2.1 mm I.D, Metal free hardware)
Column Cat. No. : 5020-76210
Eluent : A) CH₃CN
 B) 200 mM HCOONH₄ in H₂O

Time (min)	A (Vol. %)	B (Vol. %)
0	70	30
30	45	55

Analyte

1. TATGACGTTCTGATGCT 18 mer
 2. CATGACGTTCTGATGCT 18 mer
 3. TCTATGACGTTCTGATGCT 20 mer
 4. GCTATGACGTTCTGATGCT 20 mer
- (20 mmol/mL of each)

Flow Rate : 0.2 mL/min
Col. Temp. : 40 $^{\circ}$ C
Detection : UV 260 nm
Injection Vol. : 1 μ L