


InertSearch for GC Data No. GA014-0000

## **Application of the dioxins**

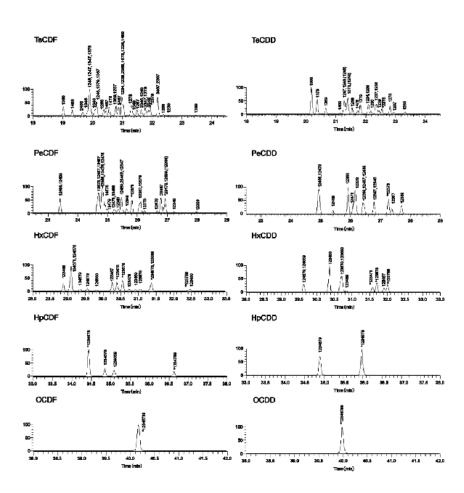

System : GC/MS

Column : InertCap 5MS/Sil

0.25mmI.D. × 60m df=0.1 $\mu$ m

Column Temp.  $:120^{\circ}\text{C}(1\text{min})\rightarrow 20^{\circ}\text{C}/\text{min}\rightarrow 200^{\circ}\text{C}(5\text{min})$ 

 $\rightarrow$ 3°C/min $\rightarrow$ 300°C(3min)

Carrier Gas : He
Detection : MS SIM
Sample : Dioxins