


InertSearch™ for GC

InertCap® Applications

Organic solvents

Data No. GA027-0000

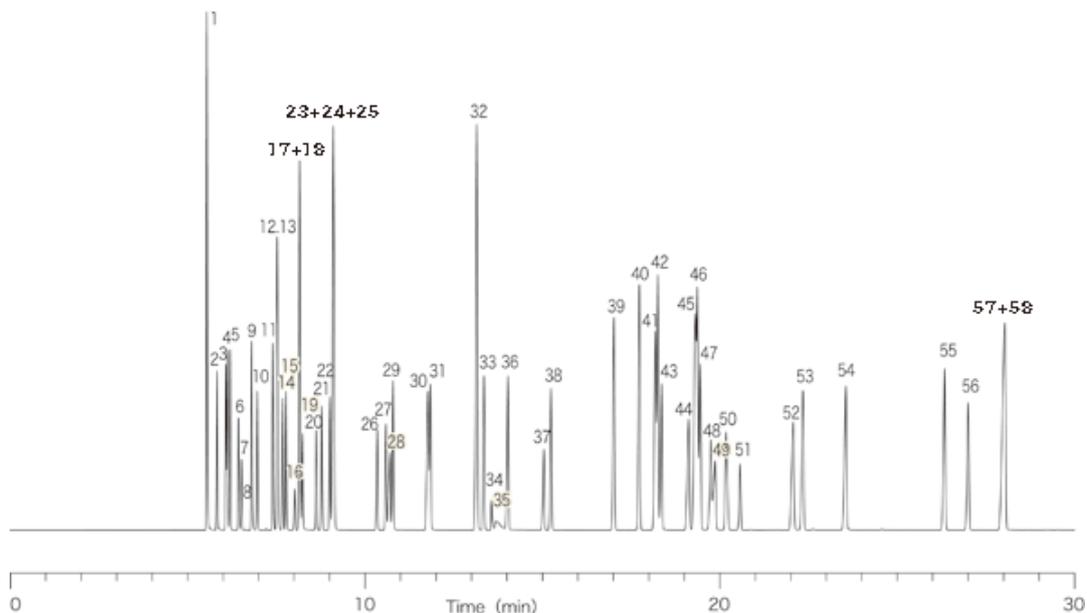

System	: GC/FID	Solutes : 1. Methanol	26. Trichloroethylene	51. 1,1,2,2-tetrachloroethane
Column	: InertCap 5MS/Sil 0.25mm I.D. × 60m df=0.25µm	2. Ethanol	27. 1,4-Dioxane	52. Methylcyclohexanol
Column Temp.	: 40°C (5min) → 4°C/min → 230°C	3. Acetone	28. Ethyl cellosolve	53. Methylcyclohexanone
Carrier Gas	: He	4. i-Propanol	29. n-Propyl acetate	54. Phenol
Injection	: Split	5. Ethyl ether	30. i-Amyl alcohol	55. o-Dichlorobenzene
Detection	: FID	6. Methyl acetate	31. Methyl i-butyl ketone	56. o-Cresol
Sample	: Organic solvents	7. Dichloromethane	32. n-Amyl alcohol	57. m-Cresol
		8. Carbon disulfide	33. Toluene	58. p-Cresol
		9. n-Propanol	34. i-Butyl acetate	
		10. trans-1,2-Dichloroethylene	35. N,N-Dimethylformamide	
		11. Methyl ethyl ketone	36. Methyl n-butyl ketone	
		12. Hexane	37. n-Butyl acetate	
		13. 2-Butanol	38. Tetrachloroethylene	
		14. cis-1,2-Dichloroethylene	39. Chlorobenzene	
		15. Ethyl acetate	40. Ethylbenzene	
		16. Chloroform	41. m-Xylene	
		17. i-Butanol	42. p-Xylene	
		18. Tetrahydrofuran	43. i-Amyl acetate	
		19. Methyl cellosolve	44. Cyclohexanol	
		20. 1,1,1-Trichloroethane	45. Styrene	
		21. 1,2-Dichloroethane	46. Cyclohexanone	
		22. i-Propyl acetate	47. o-Xylene	
		23. Benzene	48. Cellosolve acetate	
		24. n-Butanol	49. Butyl cellosolve	
		25. Carbon Tetrachloride	50. n-Amyl acetate	