


ϵ -Caprolactam

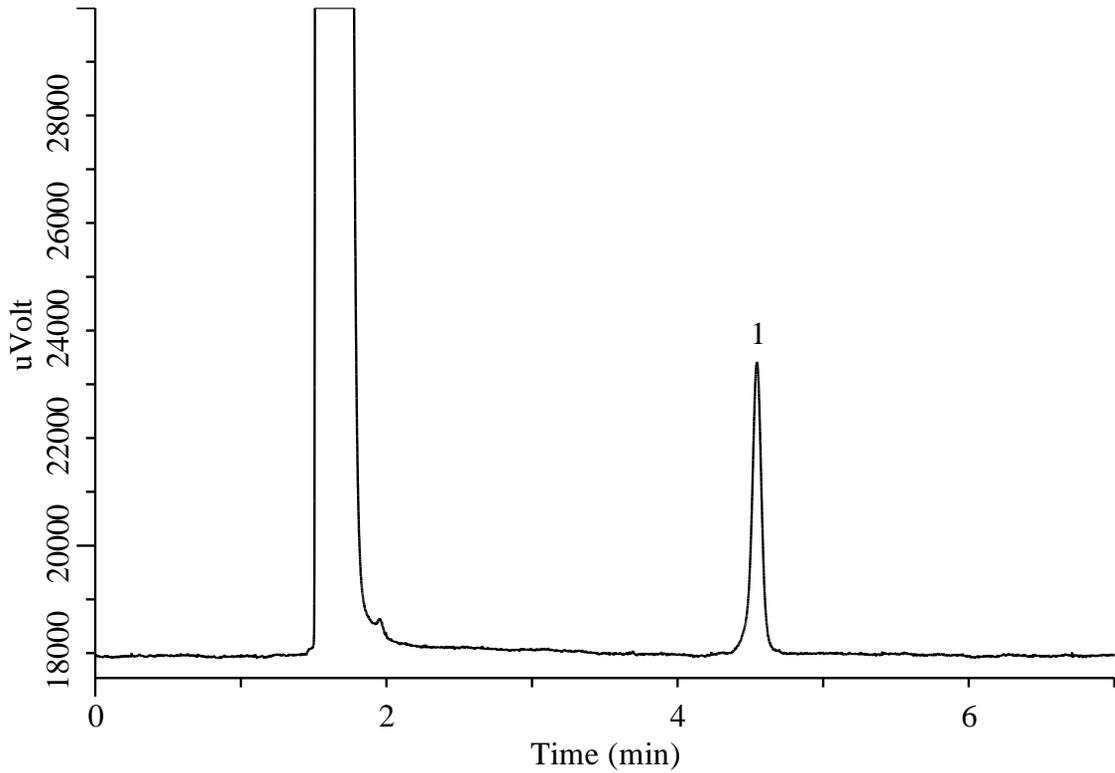

Conditions

System	: GC/FID
Column	: InertCap 1 0.32 mm I.D. x 30 m df = 5.00 μ m
Col.Cat. No.	: 1010-11249
Col. Temp.	: 240 $^{\circ}$ C
Carrier Gas	: He 110 kPa
Injection	: Split flow 18 mL/min 240 $^{\circ}$ C
Detection	: FID Range 10 ⁰ 240 $^{\circ}$ C
Sample Size	: 15 mg/L in 20% Ethanol solution 1.0 μ L
Analyte	: 1. ϵ -Caprolactam