


PTH-Amino acids

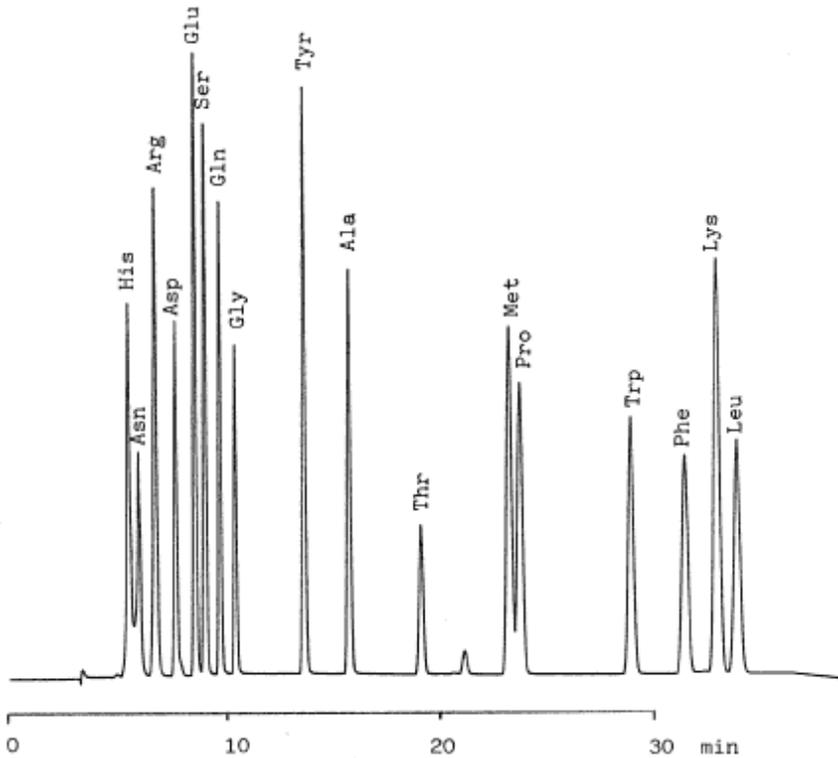

Column : Inertsil ODS-2
5µm
250 × 4.6mmI.D.

Column Cat. No. : 5020-01102

Eluent : A: CH₃CN
B: H₂O

A/B = 20/80 - 5 min - 30/70 - 25 min - 40/60

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254nm

Sample : PTH-amino acids

Solutes

- : 1. Histidine (His)
2. Asparagine (Asn)
3. Arginine (Arg)
4. Asparatic acid (Asp)
5. Glutamic acid (Glu)
6. Serine (Ser)
7. Glutamine (Gln)
8. Glycine (Gly)
9. Tyrosine (Tyr)
10. Alanine (Ala)
11. Threonine (Thr)
12. Methionine (Met)
13. Proline (Pro)
14. Tryptophan (Trp)
15. Phenylalanine (Phe)
16. Lysine (Lys)
17. Leucine (Leu)