


Free Amino acids and Organic acids

Data No. LA010-0000

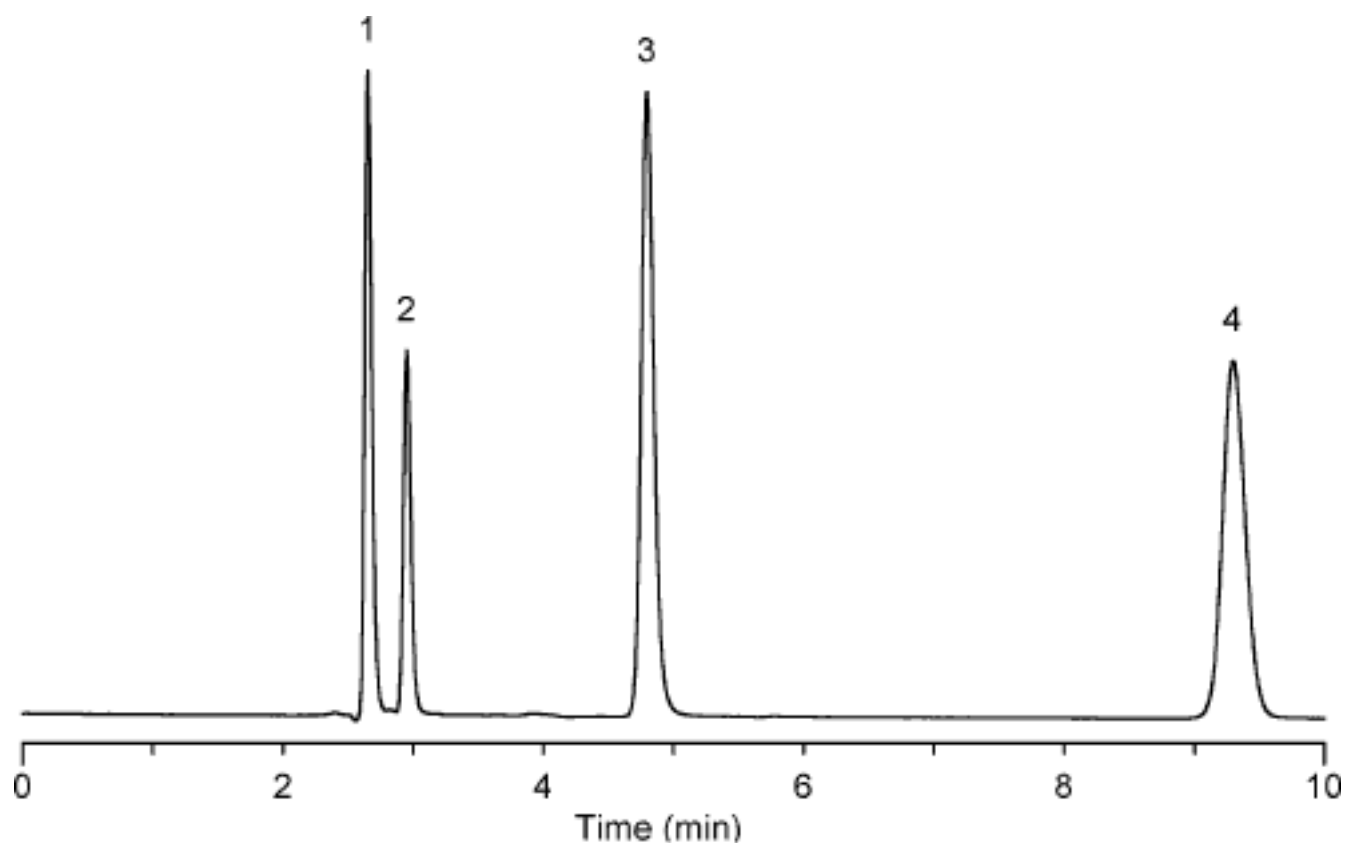

Column : Inertsil CN-3
5 μ m
250 \times 4.6 mm I.D.

Column Cat. No. :5020-01941

Eluent : 20 mM KH₂PO₄ (pH 4.0; H₃PO₄)

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 5 μ L

Sample : Free amino acids

Solutes : 1. Asparagine \cdot H₂O (0.75 mg/mL)
2. Aspartic acid (0.75 mg/mL)
3. Fumaric acid (0.01 mg/mL)
4. Maleic acid (0.01 mg/mL)