


## Vitamins

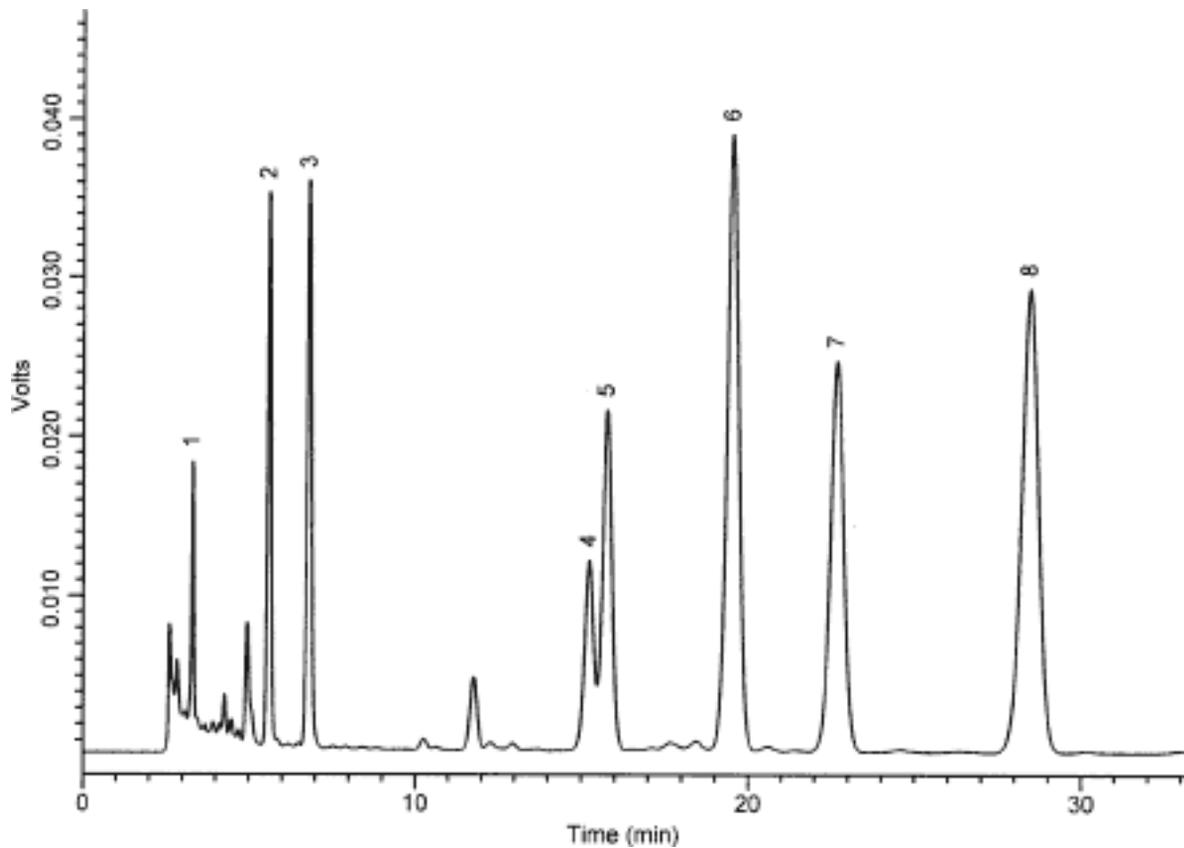

Column	:Inertsil ODS-3 250 × 4.6 mm I.D.	Solutes	: 1. Vitamin K3 (37 µg/mL) 2. Vitamin A (74 µg/mL) 3. Vitamin A acetate (222 µg/mL) 4. Vitamin D2 (111 µg/mL) 5. Vitamin D3 (37 µg/mL) 6. $\alpha$ -Tocopherol (667 µg/mL) 7. Vitamin E acetate (667 µg/mL) 8. Vitamin K1 (370 µg/mL)
Column Cat. No. :			
Eluent	:CH <sub>3</sub> CN/CH <sub>3</sub> OH = 75/25		
Flow rate	:1.0 mL/min.		
Col. Temp.	:40°C		
Detection	:UV 280 nm		
Injection Vol.	:5 µL		
Sample	:Fat-soluble Vitamins		