


InertSearch™ for LC

Inertsil® Applications

Carotene

Data No. LA026-0000

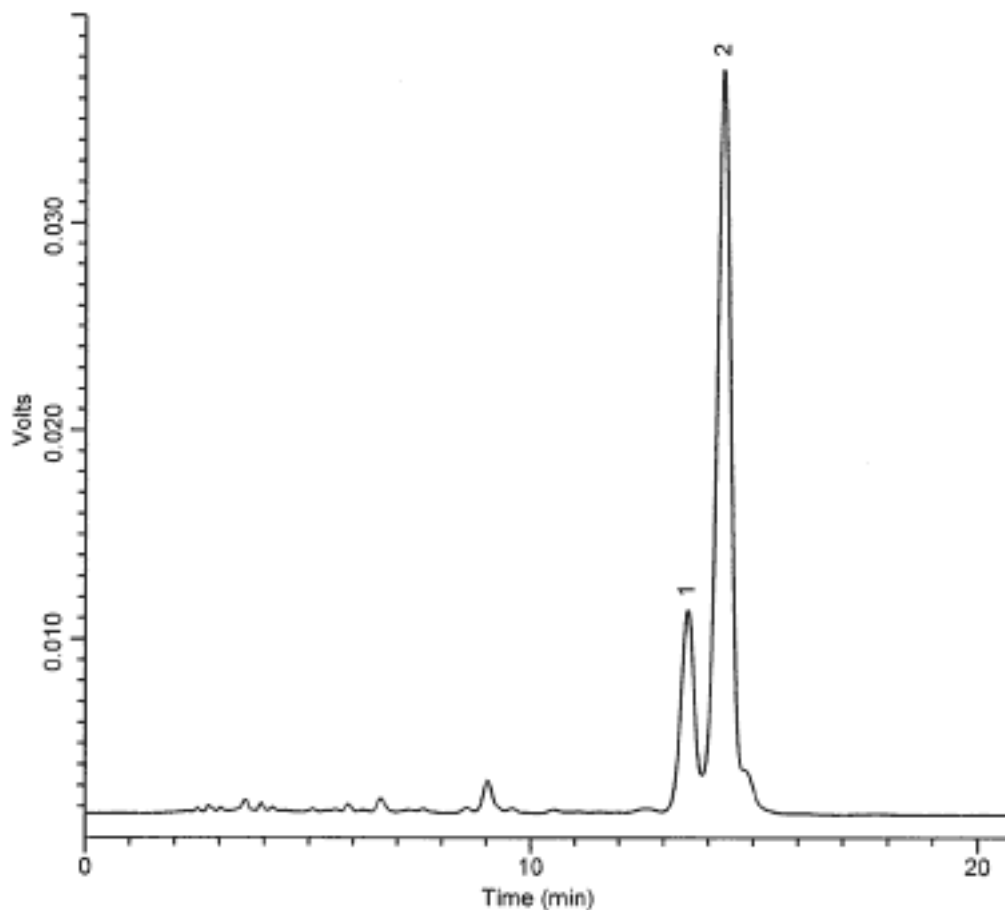

Column : Inertsil ODS-3
150 × 4.6 mm I.D.

Solutes : 1. α -Carotene (0.022 mg/mL)
2. β -Carotene (1.935 mg/mL)

Column Cat. No. :

Eluent : CH₃CN/Ethanol = 70/30

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 450 nm

Injection Vol. : 2 μ L

Sample : Carotenes