


InertSearch™ for LC

Inertsil® Applications

Vitamins

Data No. LA033-0000

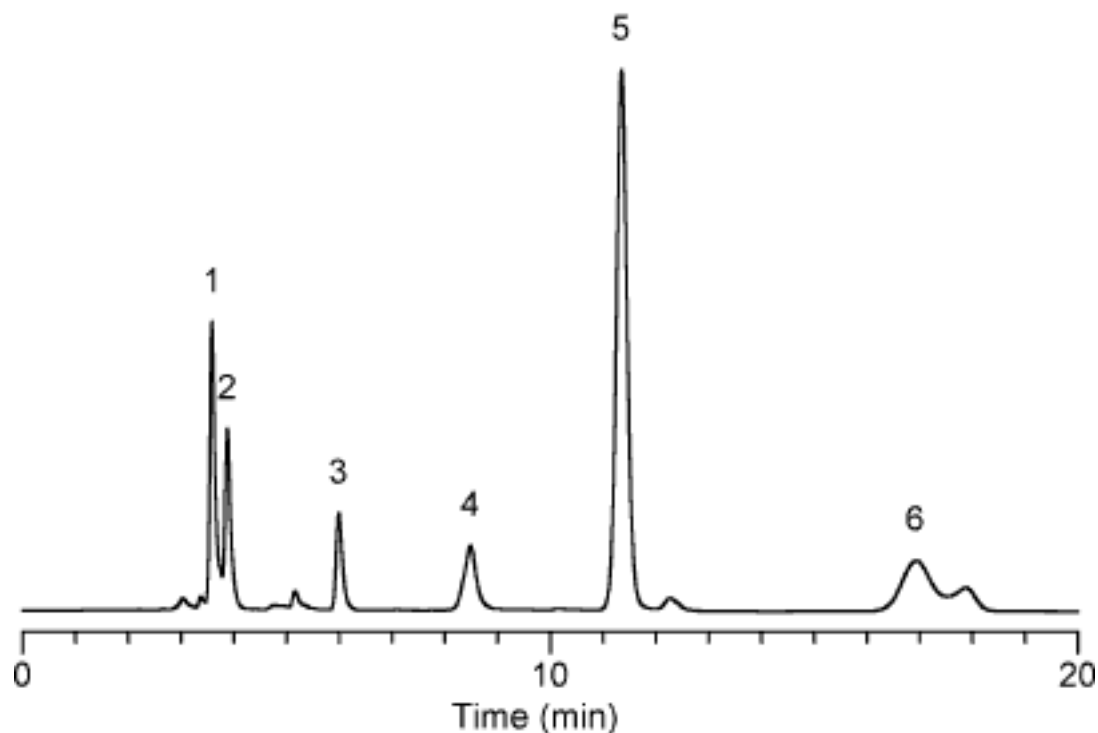

Column	: Inertsil CN-3 5µm 250 × 4.6 mm I.D.	Solutes	: 1. Vitamin K1 (0.14 mg/mL, tR 3.587) 2. Vitamin A acetate (0.14 mg/mL, tR 3.880) 3. Vitamin K3 (0.14 mg/mL, tR 5.993) 4. Vitamin E (0.14 mg/mL, tR 8.487) 5. Vitamin D3 (0.14 mg/mL, tR 11.347) 6. Vitamin A (0.14 mg/mL, tR 16.920)
Column Cat. No. :			
Eluent	: n-Hexane/Ethanol = 99.5/0.5		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 280 nm		
Injection Vol.	: 1 µL		
Sample	: Fat-soluble Vitamins		