


Vitamins

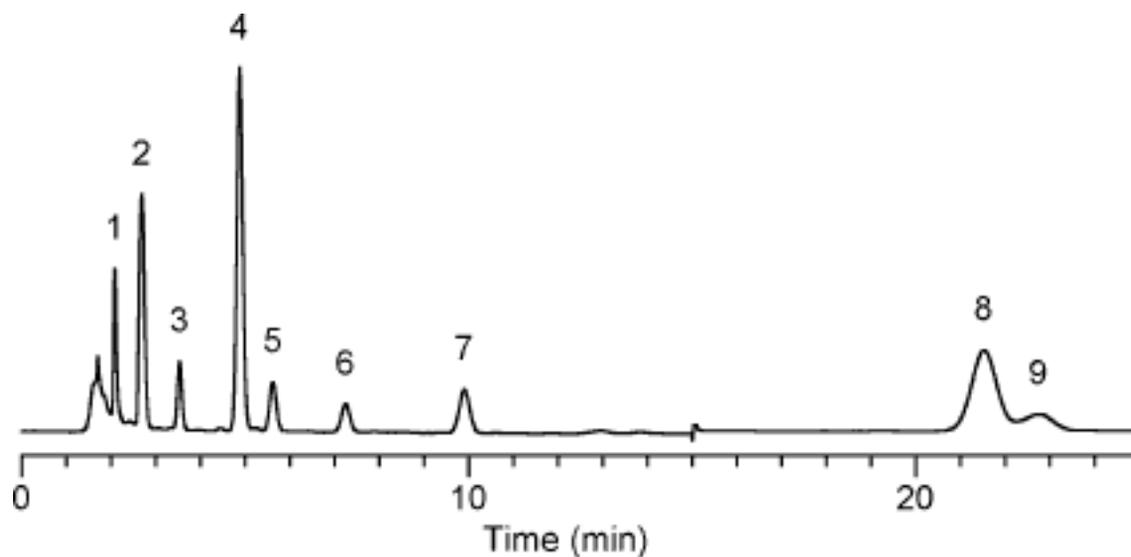

Column	: Inertsil ODS-3V 5 μ m 150 \times 4.6 mm I.D.	Solutes	: 1. Vitamin K3 (0.05 mg/mL) 2. Vitamin A (0.05 mg/mL) 3. Vitamin A acetate (0.05 mg/mL) 4. Vitamin D3 (0.05 mg/mL) 5. Vitamin E (0.05 mg/mL) 6. Vitamin E acetate (0.05 mg/mL) 7. Vitamin K1 (0.05 mg/mL) 8. α -Carotene (0.03 mg/mL) 9. β -Carotene (0.08 mg/mL)
Column Cat. No.:			
Eluent	: CH ₃ OH/Ethanol = 75/25		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: 0 - 15 min: UV 280 nm, 15 - 25 min: VIS 450 nm		
Injection Vol.	: 3 μ L		
Sample	: Fat-soluble Vitamins		