
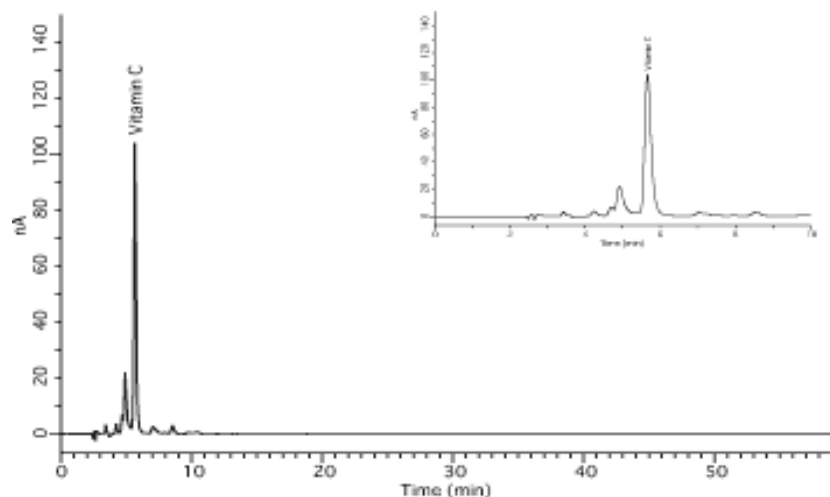


InertSearch™ for LC

Inertsil® Applications

Vitamin C

Data No. LA040-0000

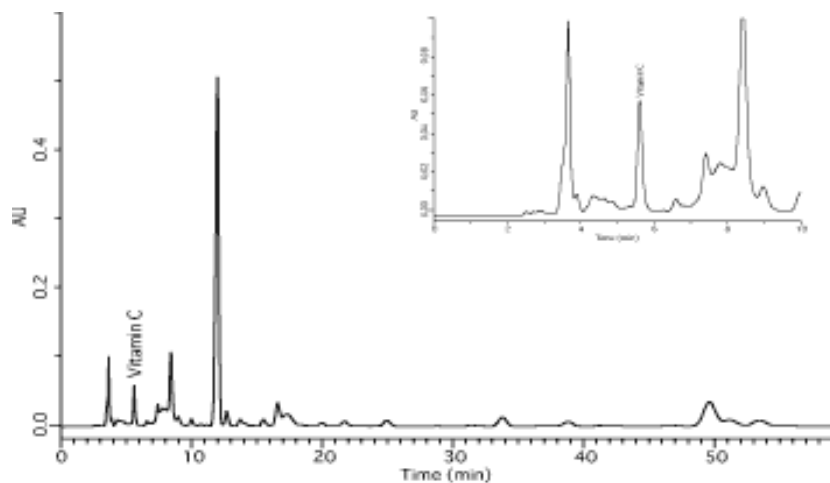

Column :Inertsil ODS-3
5 μ m
150 \times 4.6 mm I.D.

Column Cat. No. :

Eluent :100 mM KH₂PO₄, 1 mM EDTA (pH 3.0)

Flow rate :0.6 mL/min

Col. Temp. :25°C

Detection :UV 254 nm
ECD 200 mV vs. Ag/AgCl

Injection Vol. :10 μ L

Sample :Liquid-liquid extract of urine

Solutes :Vitamin C