


InertSearch for LC Data No. LA054-0000

## Analysis of Vitamin D2 and D3

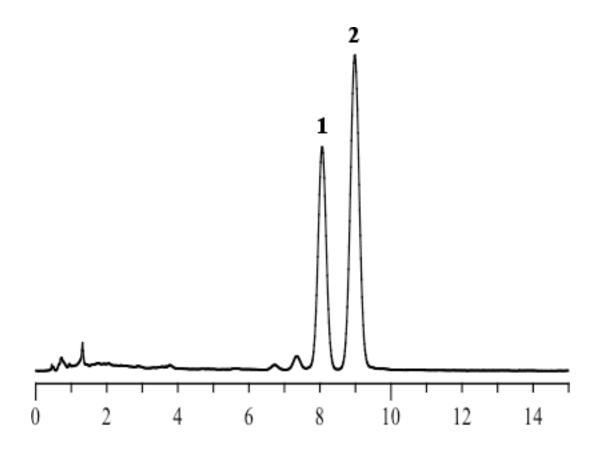

Solutes

: 1. Vitamin D2

2. Vitamin D3

Column : Inertsil ODS-P

 $5 \, \mu m$ 

150 x 4.6 mm I.D.

Column Cat. No.:

Eluent : A) CH<sub>3</sub>CN

B) CH<sub>3</sub>OH

A/B = 75/25, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 280 nm

Injection Vol.

Sample :