


7. Sugars

Data No. LA062-0000

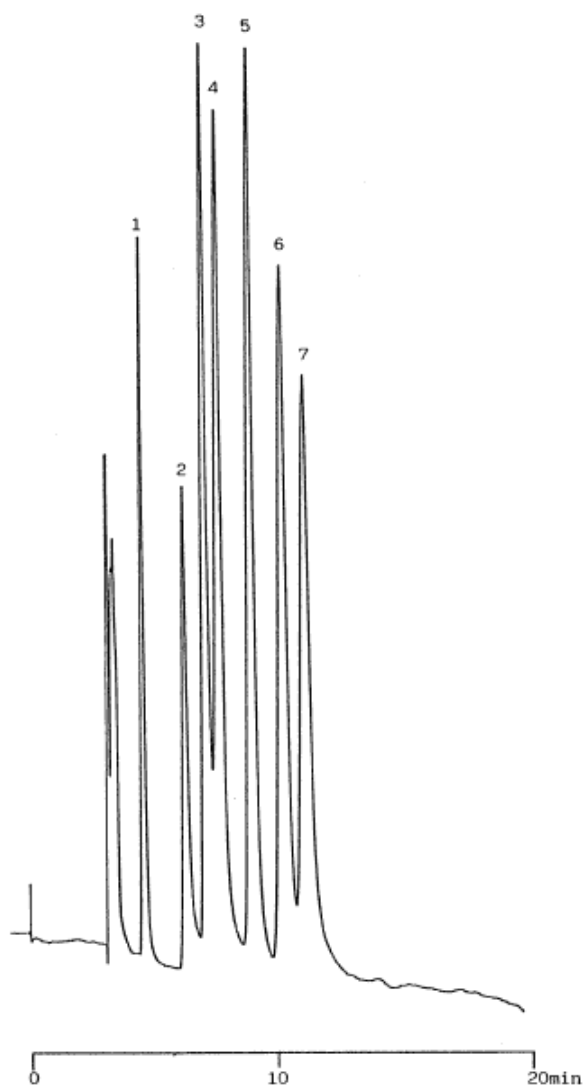

Column : Unisil Q NH2
5 μ m
250 $\times$ 4.6 mm I.D.

Column Cat. No. : 5020-03402

Eluent : CH₃CN/H₂O = 65/35

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : RI

Injection Vol. :

Sample : Sugars

Solutes : 1. Ethyleneglycol
2. Xylose
3. D (-) Fructose
4. D (+) Glucose
5. Sucrose
6. D (+) Maltose monohydrate
7. Lactose monohydrate