


Analysis of Monosaccharides and Disaccharides Data No. LA100-0000

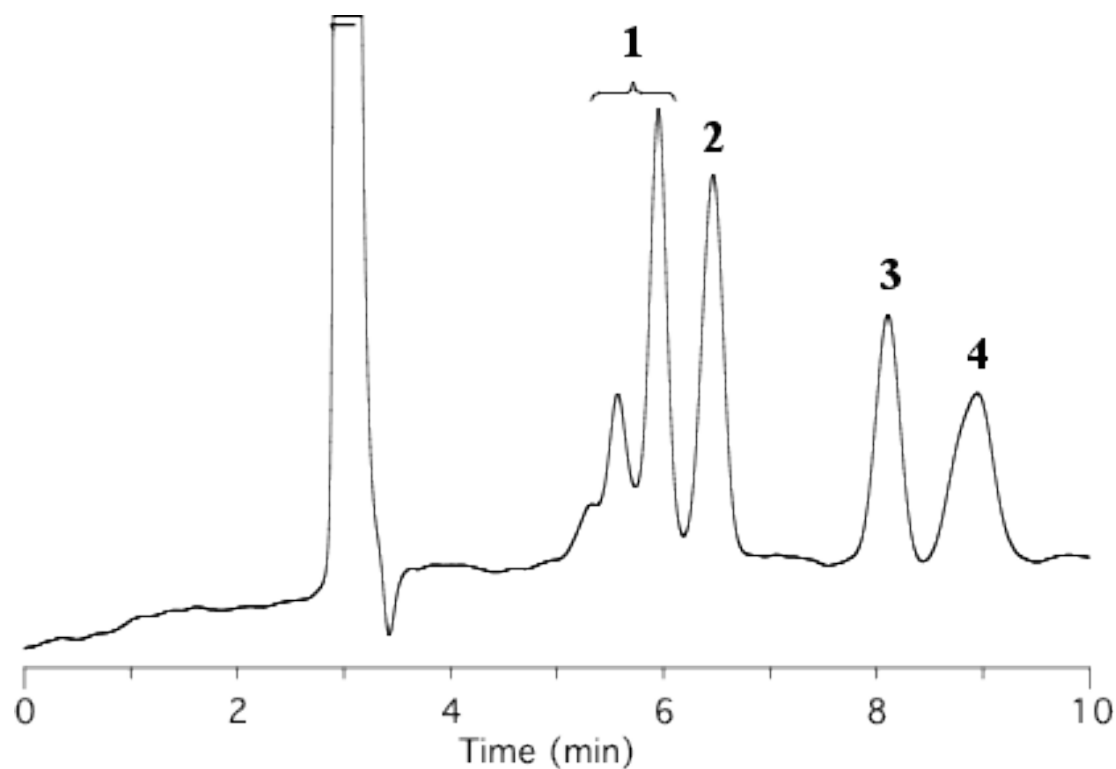

Column	: Inertsil Diol 5 µm 250 x 4.6 mm I.D.	Solutes	: 1. Fructose 2. Glucose 3. Sucrose 4. Maltose
Column Cat. No.	: 5020-05646		
Eluent	: A) CH ₃ CN B) H ₂ O A/B = 80/20, w/w		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: RI		
Injection Vol.	:		
Sample	:		