


Organic acids

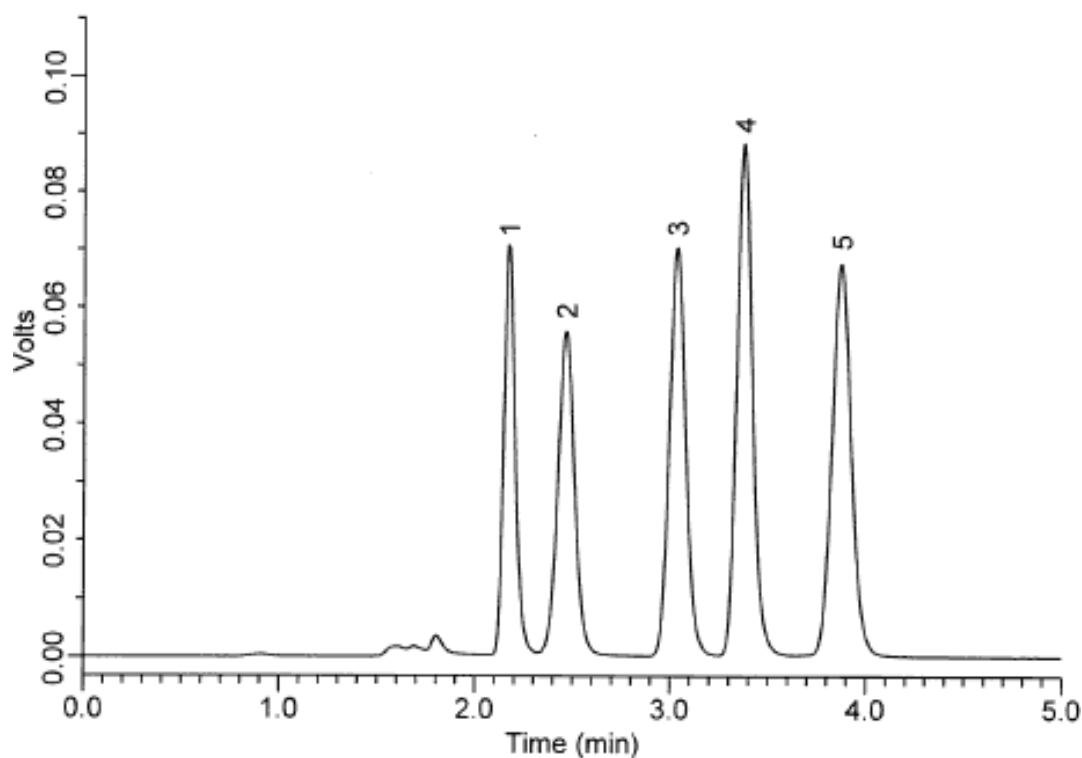

Column : Inertsil C8-3
150 × 4.6mmI.D.

Column Cat. No.: 5020-01900

Eluent : (NH₄)H₂PO₄(pH 2.8; H₃PO₄)

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. :

Sample : Organic acids

Solutes : 1. Formic acid(0.5mg/mL)
2. alpha-Ketoglutaric acid(0.05mg/mL)
3. Citric acid(0.5mg/mL)
4. Maleic acid(0.005mg/mL)
5. Fumaric acid(0.005mg/mL)