


Anion Surfactants

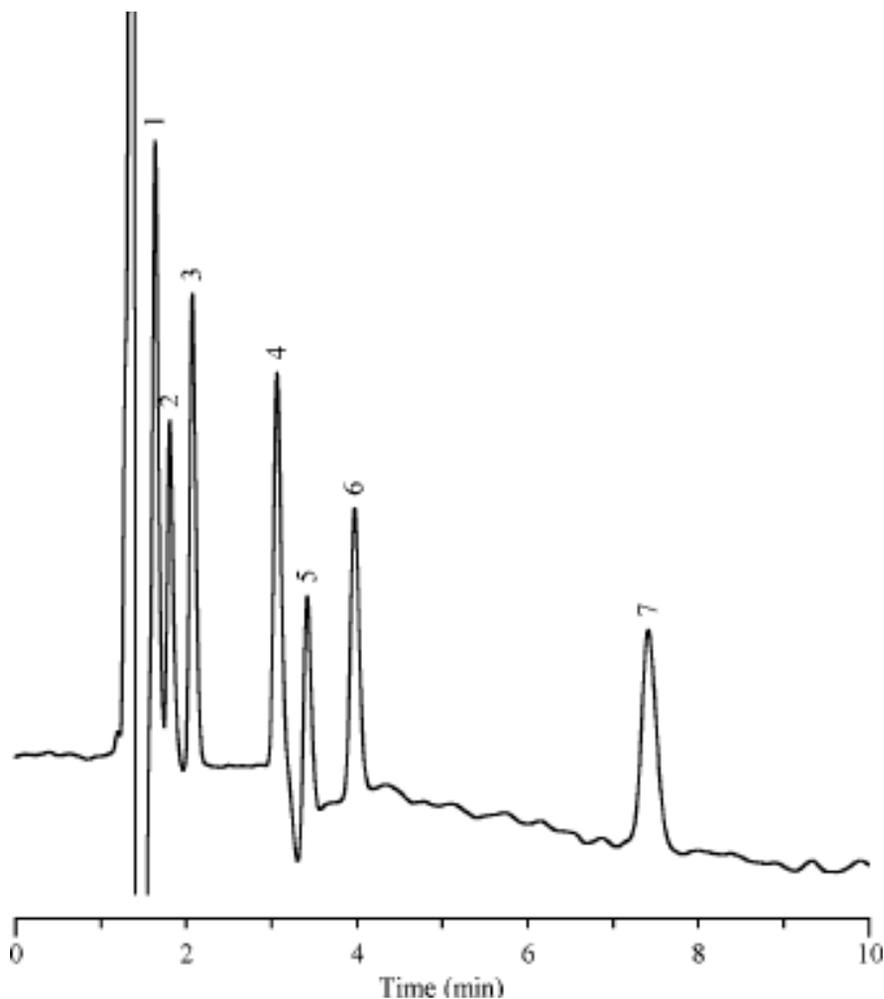

Column : Inertsil ODS-3V
 5µm
 150×4.6mmI.D.

Column Cat. No. : 5020-01801

Eluent : 0.1M NaClO₄ / acetonitrile (40/60)

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : RI

Injection Vol. :

Sample : Anion Surfactants (80mg/mL)

Solutes :1. Nonyl sulfonate Na
 2. Decyl sulfate Na
 3. Undecyl sulfonate Na
 4. Tridecyl sulfonate Na
 5. Dodecyl sulfate Sodium
 7. Tetradecyl sulfonate Na
 8. Hexadecyl sulfonate Na