


Tartaric acid

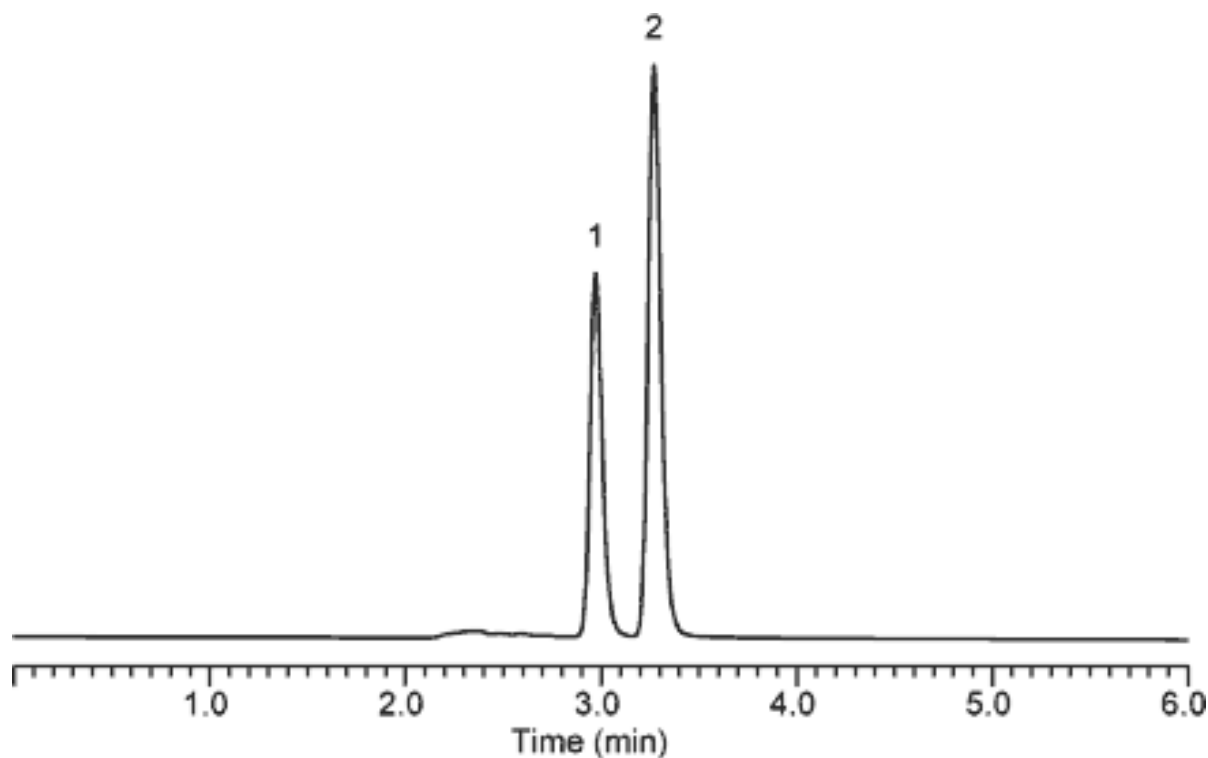

Column	: Inertsil ODS-3V 5 μ m 250 $\times$ 4.6mm I.D.	Solutes	: 1. meso-Tartaric acid 2. DL-Tartaric acid
Column Cat. No.	: 5020-01802		
Eluent	: 50mM Phosphate buffer (pH2.5)		
Flow rate	: 1mL/min		
Col. Temp.	: 40°C		
Detection	: UV 210 nm		
Injection Vol.	: 5 μ L		
Sample	: meso-Tartaric acid 1-hydrate (0.5 mg/mL) DL-Tartaric acid (0.5 mg/mL)		