


## InertSearch for LC Data No. LA127-0000

## **Organic acids**

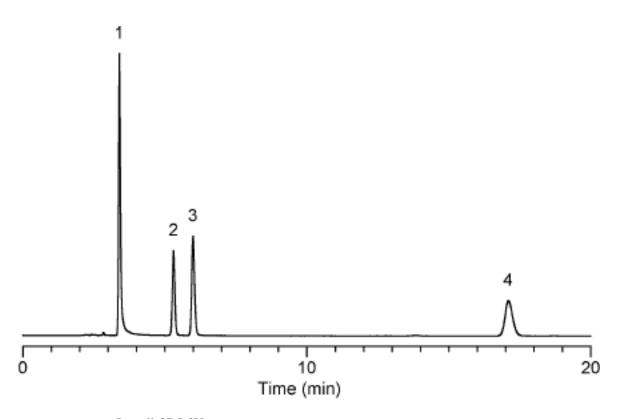

Column : Inertsil ODS-3V

5µm

250×4.6mmI.D.

Column Cat. No.: 5020-01802

Eluent : 20 mM KH2PO4 in H2O (pH 3.0, H3PO4)

Flow rate : 1.0mL/min

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 20µL

Sample : Organic acids

Solutes :1. Formic acid (0.1 %)

2. Lactic acid (0.1 %)

3. Acetic acid (0.1 %)

4. Propionic acid (0.1 %)