


Rhei rhizoma

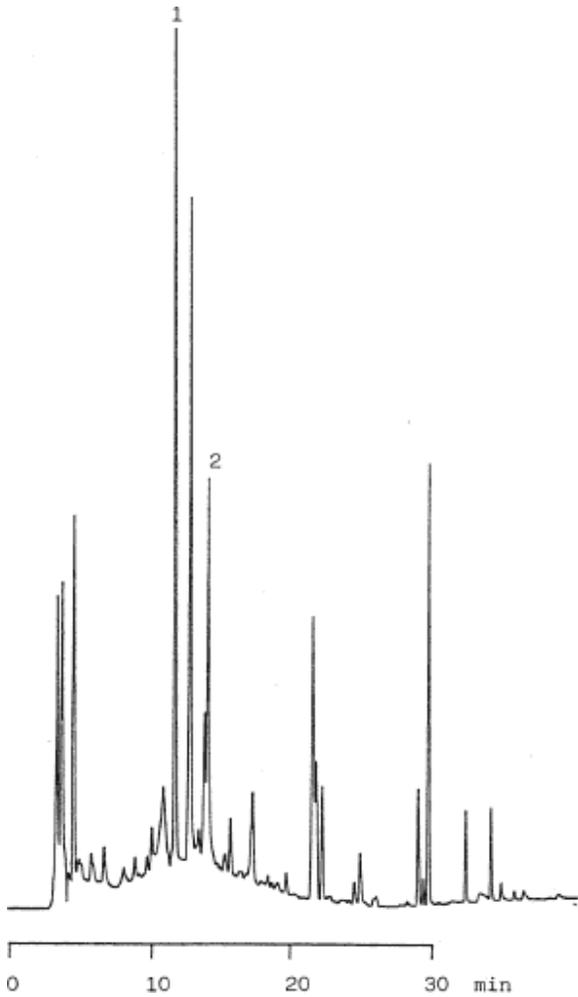

Column	: Inertsil C8 5µm 150 × 4.6mm I.D.	Solutes	: 1. Sennoside B 2. Sennoside A
Column Cat. No.	: 5020-01201		
Eluent	: A: CH ₃ CN B: 40mM NaClO ₄ A/B 15/85?(20min) → 40/60 → (10min)?90/10		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		