


InertSearch™ for LC

Inertsil® Applications

Zingiberis rhizoma

Data No. LA132-0000

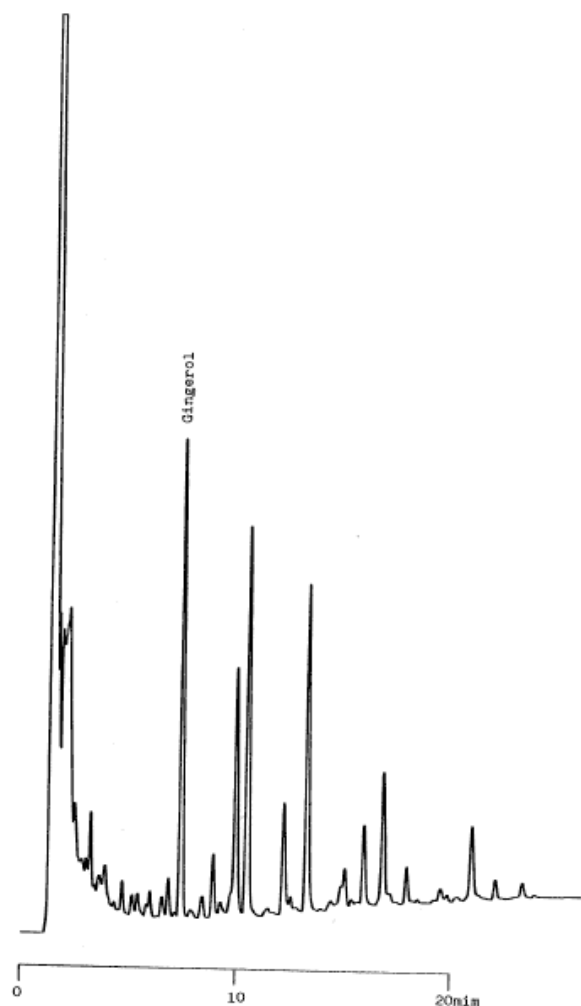

Column : Inertsil ODS-2
5 μ m
150 \times 4.6 mm I.D.
Column Cat. No. : 5020-01101
Eluent : A:CH₃CN
B:H₂O
A/B 50/50 \rightarrow (20min) \rightarrow 90/10
Flow rate : 1.0 mL/min.
Col. Temp. : 40°C
Detection : UV 254 nm
Injection Vol. :
Sample :
Solutes : Gingerol