


Herb medicine

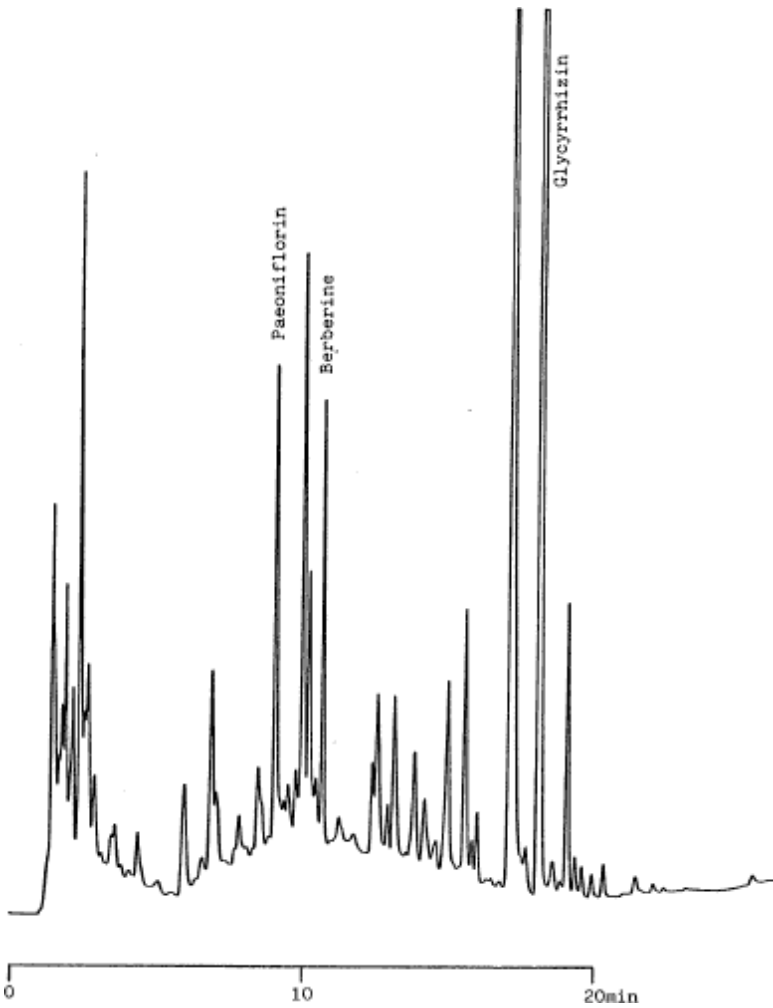

Column : Inertsil ODS-2
 5 μ m
 150 \times 4.6 mmI.D.
 Column Cat. No.: 5020-01101
 Eluent : A:CH₃CN
 B:0.05% H₃PO₄
 A/B 10/90 \rightarrow (25min) \rightarrow 90/10
 Flow rate : 1.0 mL/min.
 Col. Temp. : 40°C
 Detection : UV 254 nm
 Injection Vol. :
 Sample :
 Solutes : Paeoniflorin
 Berberin
 Glycyrrhizin