


## Puerariae radix

Data No. LA134-0000

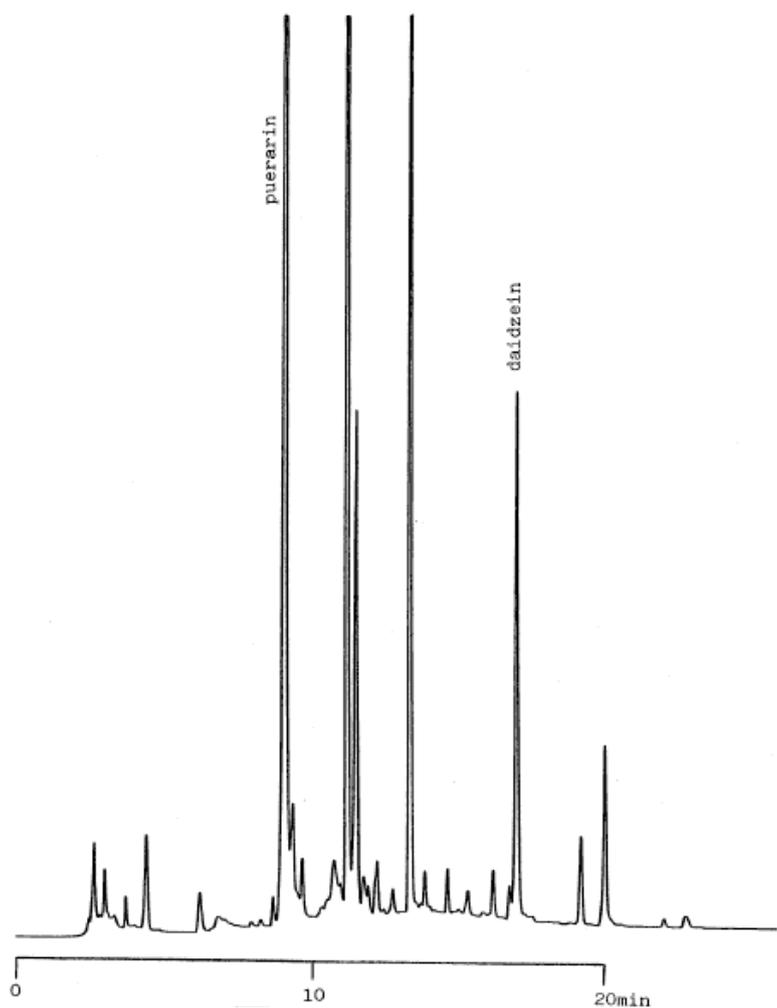

Column : Inertsil ODS-2  
5 $\mu$ m  
250  $\times$  4.6 mm I.D.  
Column Cat. No. : 5020-01102  
Eluent : A:CH<sub>3</sub>CN  
B:H<sub>2</sub>O  
A/B 10/90 $\rightarrow$ (10min) $\rightarrow$ 60/40  
Flow rate : 1.0 mL/min.  
Col. Temp. : 40°C  
Detection : UV 254 nm  
Injection Vol. :  
Sample :  
Solutes :Puerarin  
Daidzein