


InertSearch™ for LC

Inertsil® Applications

Bilirubin in Crude drug

Data No. LA153-0000

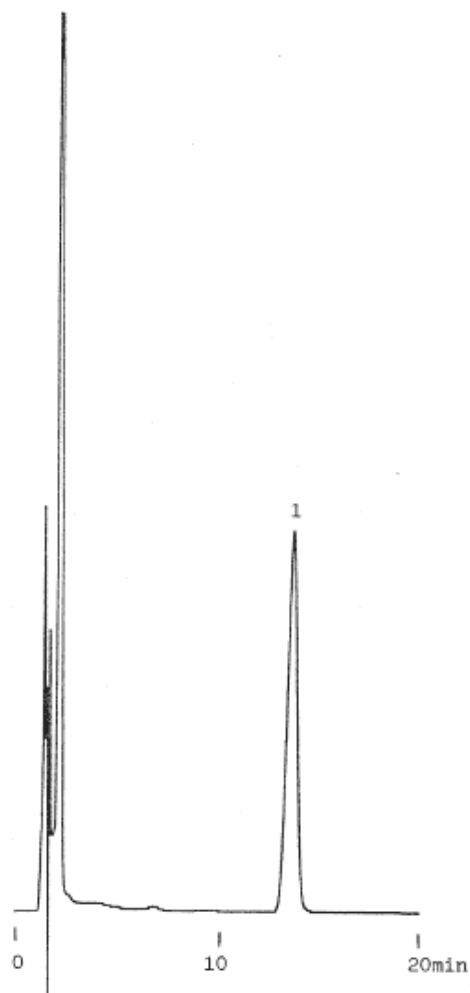

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Bilirubin

Column Cat. No. : 5020-13001

Eluent : CH₃CN/H₂O/CF₃COOH=900/100/1

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : VIS 450 nm

Injection Vol. :

Sample : 1. Bilirubin