


## InertSearch for LC Data No. LA173-0000

## **Analysis of Bisacodyl**

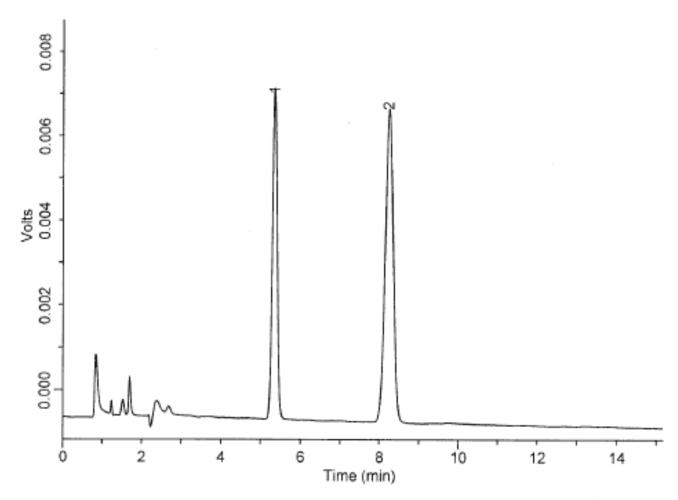

Solutes

Rs=8.8

: 1. Ethyl paraben(1.5µg/mL)

2. Bisacodyl(1.0µg/mL)

Column

:Inertsil ODS-3V

5µm

250 × 4.6mmI.D.

Column Cat. No.:

Eluent

:CH3OH/CH3CN/10mM Citric acid=1/1/2

Flow rate

:2.0 mL/min.

Col. Temp.

:30°C

Detection

:UV 254 nm

Injection Vol.

:20µL

Sample

:1. Ethyl paraben(1.5µg/mL)

2. Bisacodyl(1.0µg/mL)