


Analysis of Betamethasone sodium phosphate

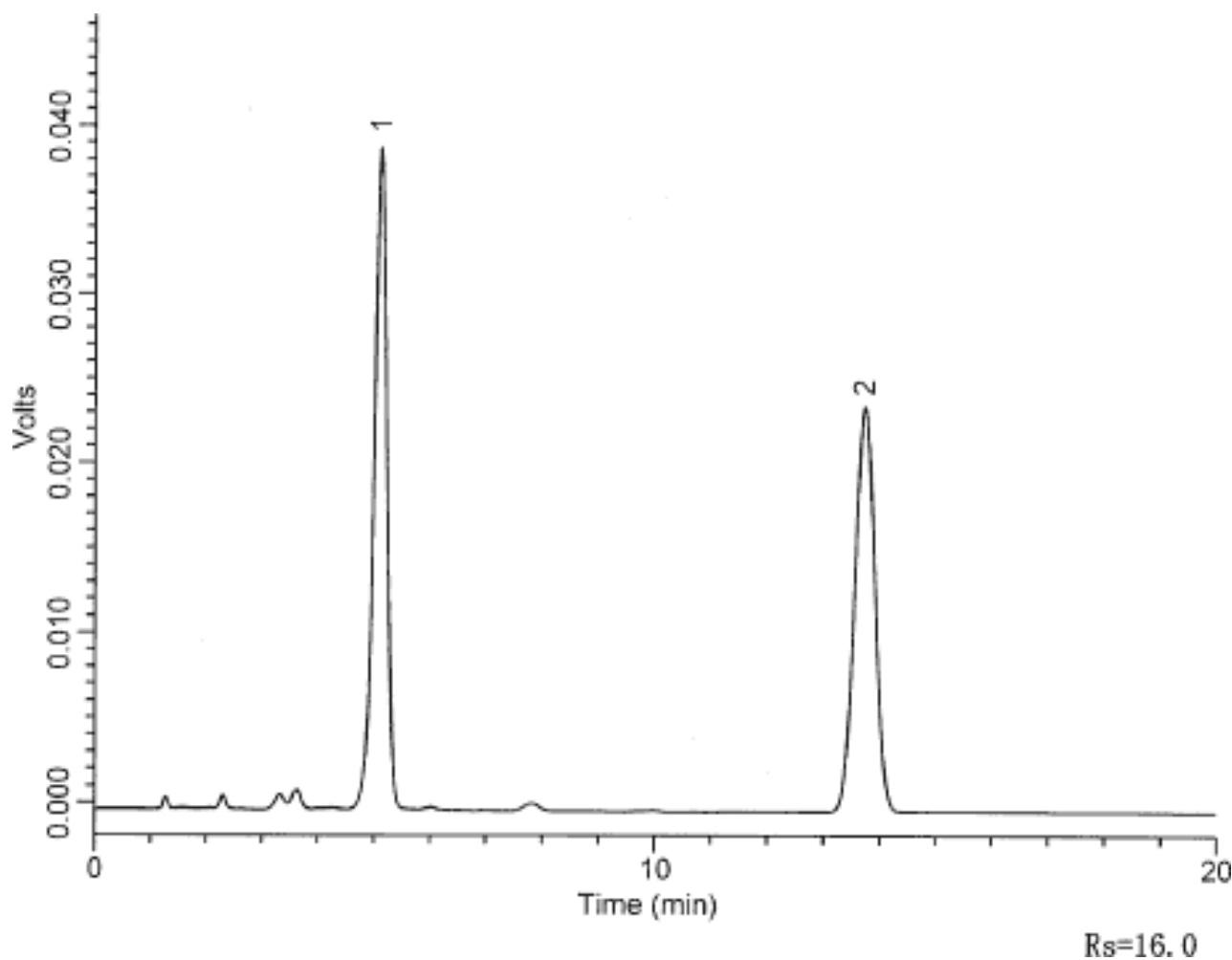

Column	: Inertsil ODS-3V 5 μ m 250 $\times$ 4.6mm I.D.	Solutes	: 1. Betamethasone sodium phosphate (0.1 mg/mL) 2. n-Butyl paraben (0.02mg/mL)
Column Cat. No. :			
Eluent	: CH ₃ OH/105mM (n-C ₄ H ₉) ₄ NBr+51mM KH ₂ PO ₄ +22.5mM Na ₂ HPO ₄ mM Na ₂ HPO ₄ =60/40		
Flow rate	: 2.0 mL/min.		
Col. Temp.	: 25°C		
Detection	: UV 254 nm		
Injection Vol.	: 10 μ L		
Sample	: 1. β -methasone sodium phosphate (0.1mg/mL) 2. n-Buthyl paraben (0.02mg/mL)		