


Analysis of Alprostadi

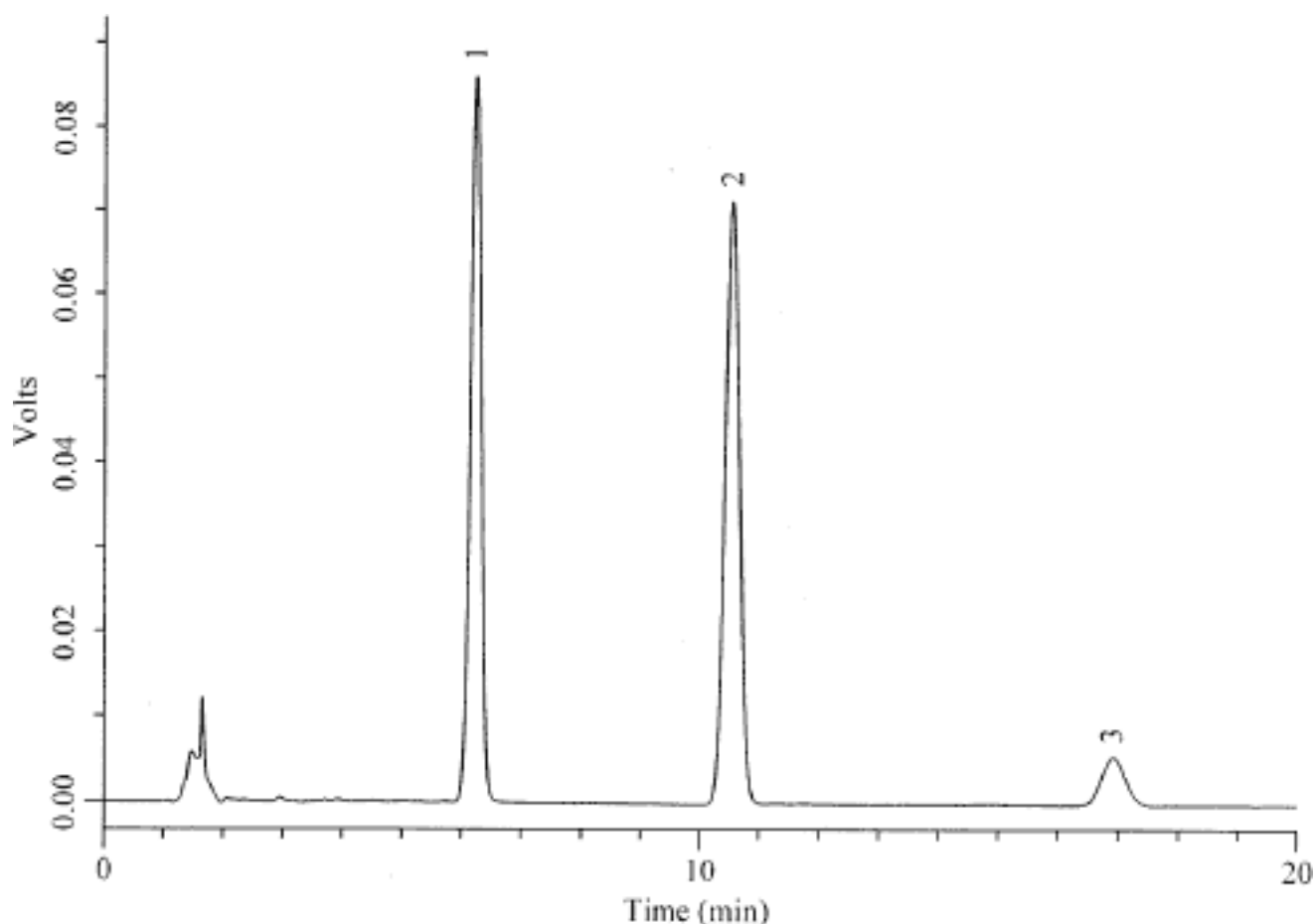

Column	:Inertsil ODS-3V 5µm 150 × 4.6mmI.D.	Solutes	: 1. Alprostadi (Prrostaglandin E1)(200µg/mL) 2. n-Propyl paraben(22µg/mL) 3. Prostaglandin A1(5µg/mL)
Column Cat. No. :			
Eluent	:CH ₃ CN/20mM KH ₂ PO ₄ =40/60		
Flow rate	:1.0 mL/min.		
Col. Temp.	:25°C		
Detection	:UV 205 nm		
Injection Vol.	:10µL		
Sample	:1. Alprostadi (Prrostaglandin E1)(200µg/mL) 2. n-Propyl paraben(22µg/mL) 3. Prostaglandin A1(5µg/mL)		