


InertSearch for LC Data No. LA193-0000

Analysis of Clofibrate

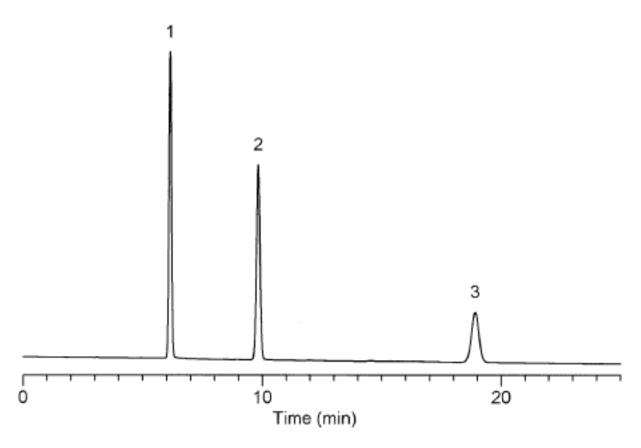

Column

:Inertsil CN-3

5µm

250 × 4.6mmI.D.

Solutes

: 1. p-Ethoxyphenol

2. p-Chlorophenol

3. Clofibrate

Column Cat. No.:

Eluent

:acetonitrile / 20mM potassium phosphate buffer

(pH3.0)=30/70

Flow rate

:1.0 mL/min.

Col. Temp.

:40°C

Detection

:UV 275 nm

Injection Vol.

:2μL

Sample

:Clofibrate (1mg/mL)

p-Chlorophenol (0.6mg/mL)

p-Ethoxyphenol (0.6mg/mL)