
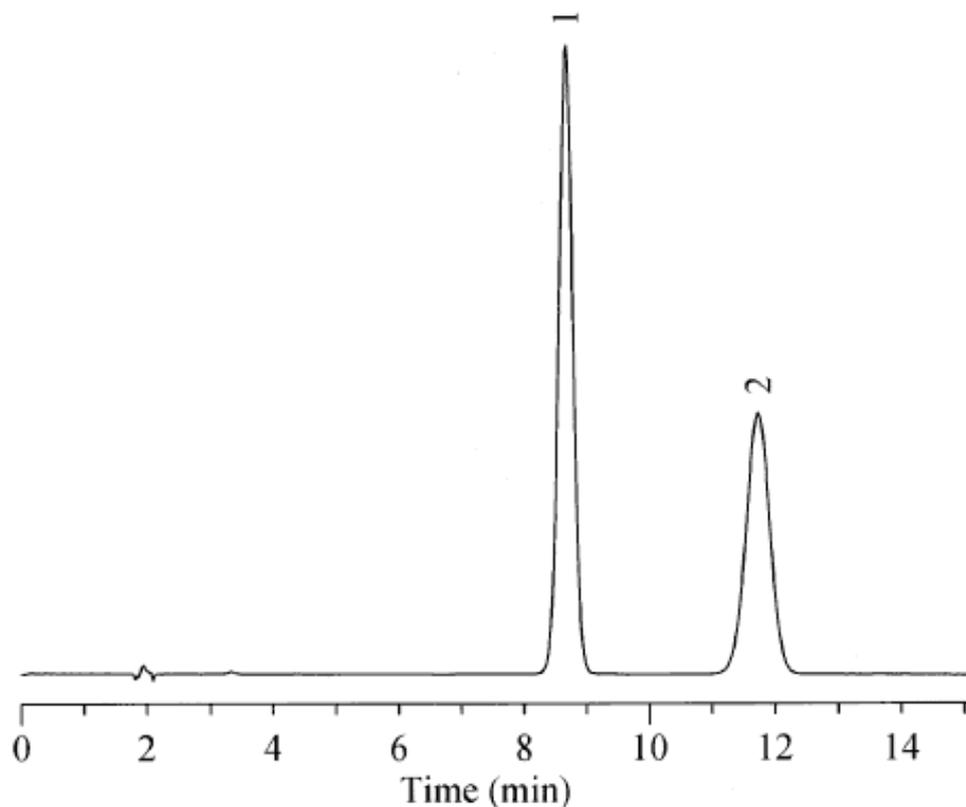


InertSearch™ for LC

Inertsil® Applications

Analysis of Carbazochrome

Data No. LA204-0000

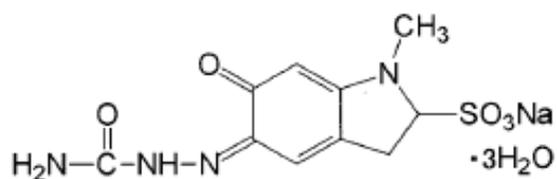

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D. Solutes : 1. Carbazochromesodium sulfonate (120ng)
2. Carbazochrome (120ng)

Column Cat. No. :

Eluent : 10mM Ammonium dihydrogen phosphate / EtOH (975 : 25)

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 360 nm

Injection Vol. :

Sample : 1. Carbazochromesodium sulfonate (120ng)
2. Carbazochrome (120ng)