


Estradiol benzoate

Data No. LA211-0000

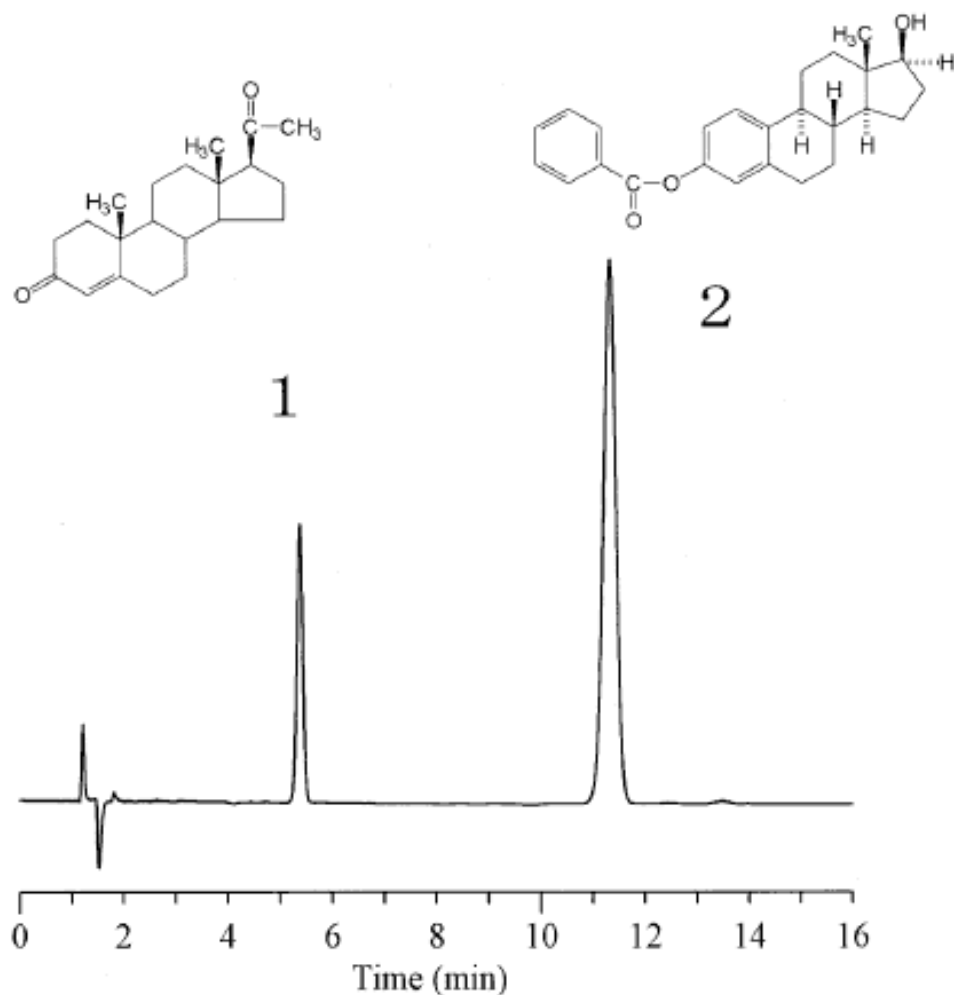

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Progesterone (0.25 μ g)
2. Estradiol benzoate (0.08 μ g)

Column Cat. No. :

Eluent : MeCN / H₂O (70 : 30)

Flow rate : 1.0 mL/min.

Col. Temp. : 35°C

Detection : UV 230 nm

Injection Vol. :

Sample : 1. Progesterone (0.25 μ g)
2. Estradiol benzoate (0.08 μ g)