


Bases of Nucleic acid and Nucleosides

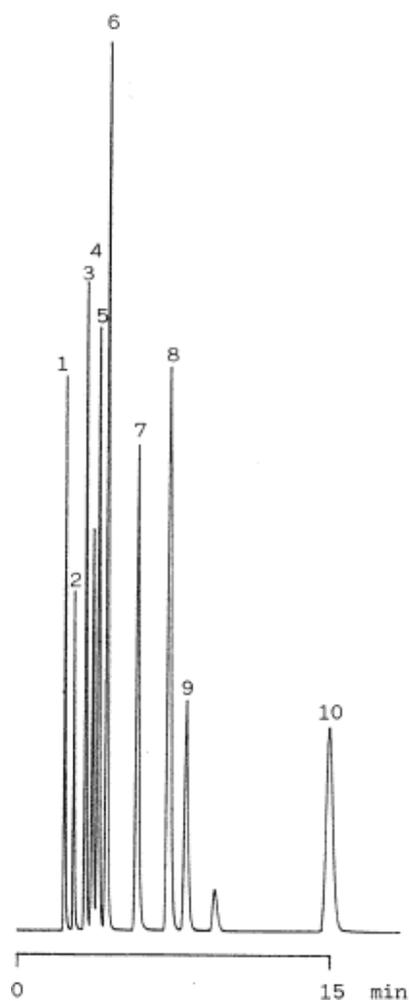

Column	:Inertsil ODS-2 5 μ m 150 \times 4.6mmI.D.	Solutes	: 1. Cytosine 2. Uracil 3. Guanine 4. Adenine 5. Cytidine 6. Uridine 7. Thymine 8. Adenosine 9. Guanosine 10. Thymidine
Column Cat. No.	:5020-01101		
Eluent	:0.1M H ₃ PO ₄ +0.2M NaC10 ₄ (pH2)		
Flow rate	:1.0 mL/min.		
Col. Temp.	:50°C		
Detection	:UV 260 nm		
Injection Vol.	:		
Sample	:		