


Analysis of gamma-Oryzanol

Data No. LA240-0000

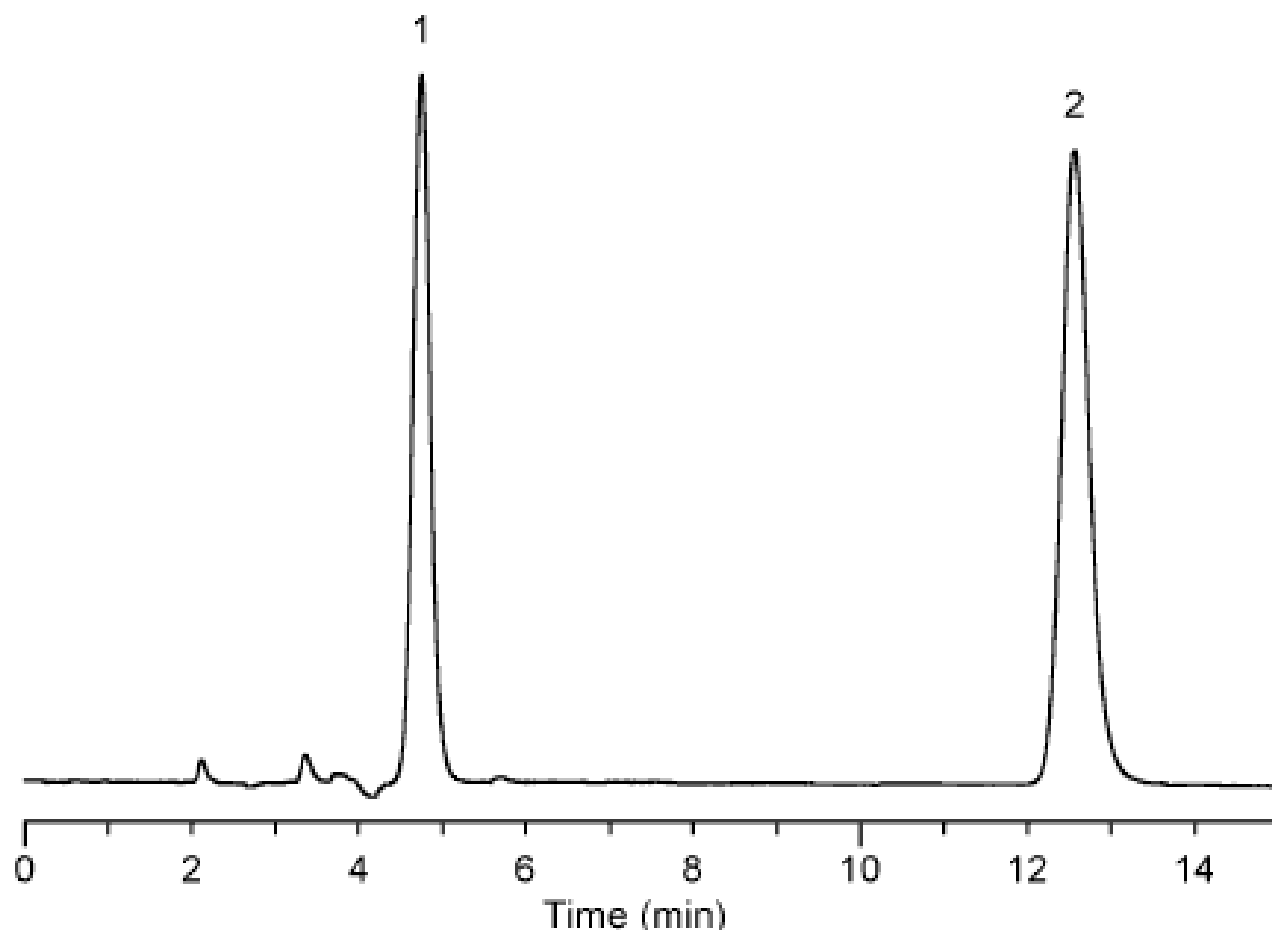

Column : Inertsil CN-3
5µm
150 × 4.6mm I.D.

Solutes : 1. gamma-oryzanol
2. vanillin

Column Cat. No. :
Eluent : hexane:2-propanol:acetic acid=94:5:1
Flow rate : 1.0mL/min
Col. Temp. : 30°C
Detection : UV 320 nm
Injection Vol. : 1µL
Sample : gamma-oryzanol (mixture 25g/mL)
vanillin (50g/mL)