


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA255-0000

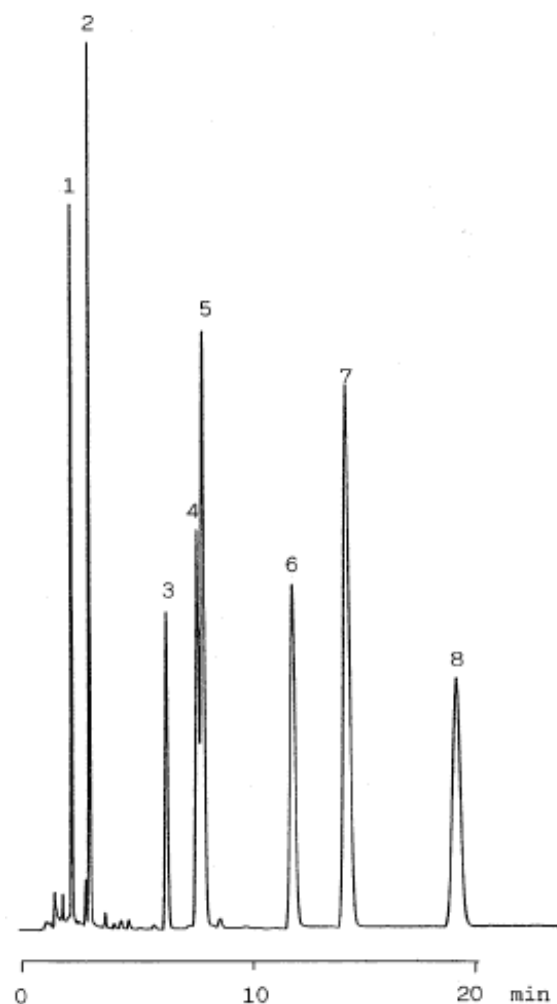

Column : Inertsil ODS-2
5 μ m
150 \times 1.5mm I.D.

Column Cat. No. :

Eluent : CH₃CN/H₂O=55/45

Flow rate : 100 μ L/min.

Col. Temp. : 50°C

Detection : UV 235nm

Injection Vol. :

Sample :

Solutes : 1. Dimethoate
2. DDVP
3. Methyl parathion
4. Malathion
5. Fenitrothion
6. Parathion
7. Diazinon
8. EPN