


# InertSearch™ for LC

Inertsil® Applications

## Siduron (cis,trans Mix)

Data No. LA260-0000

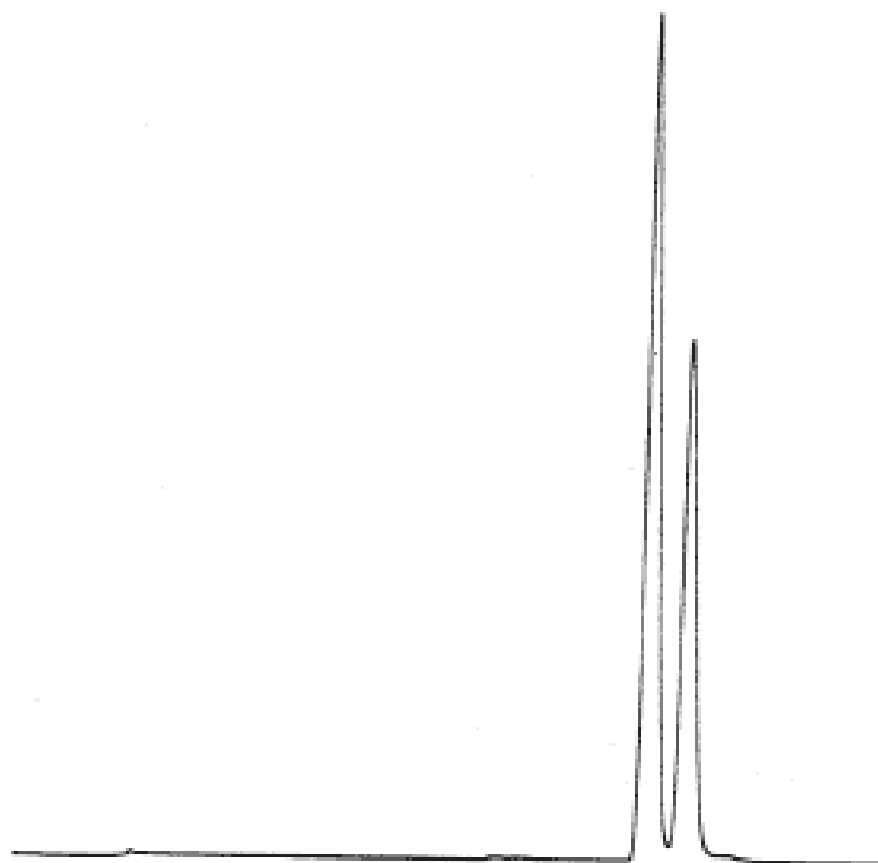

0 5 min

Column : Inertsil ODS-2  
5 $\mu$ m  
250  $\times$  4.6 mmI.D.

Column Cat. No. : 5020-01101

Eluent : CH<sub>3</sub>CN/0.1% H<sub>3</sub>PO<sub>4</sub>=50/50

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254nm

Injection Vol. :

Sample :

Solutes : Siduron