


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA269-0000

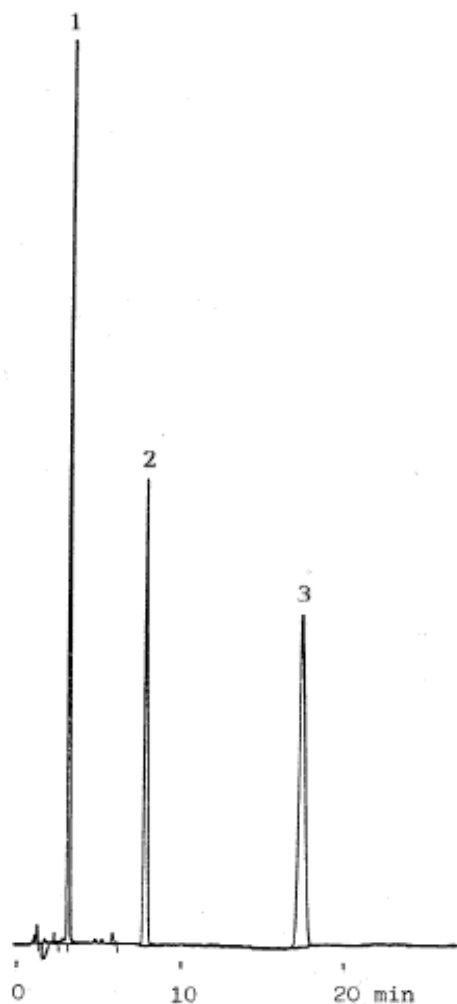

Column : Inertsil ODS-2
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. MCP
2. Pencyuron
3. Bethrodine

Column Cat. No. : 5020-01101

Eluent : CH₃CN/0.1% H₃PO₄=65/35

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :