


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA272-0000

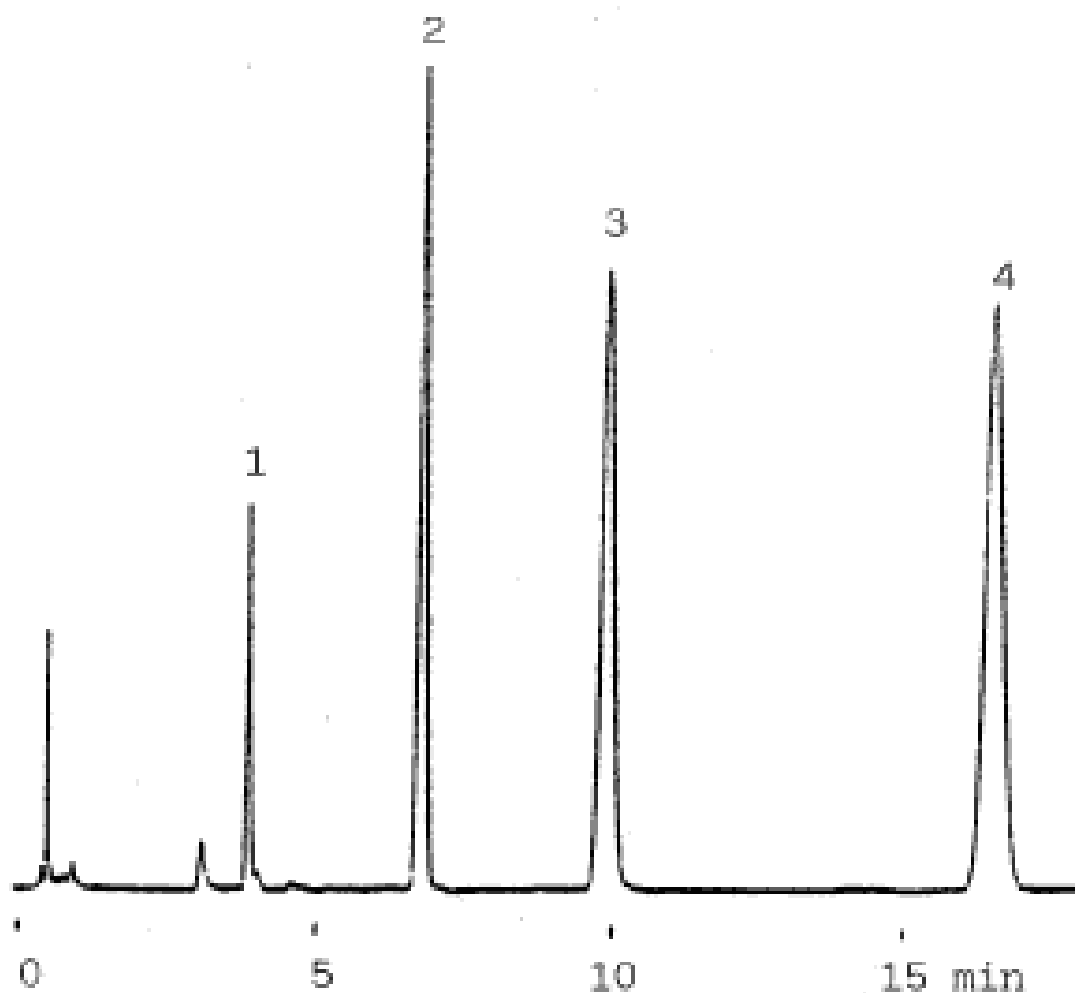

Column : Inertsil ODS-2
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/H₂O=65/35

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. :

Sample :

Solutes : 1. Fenobucarb(BPMC)
2. Bensulide
3. Isofenphos
4. Pendimethalin