


# InertSearch™ for LC

Inertsil® Applications

## Pesticides

Data No. LA279-0000

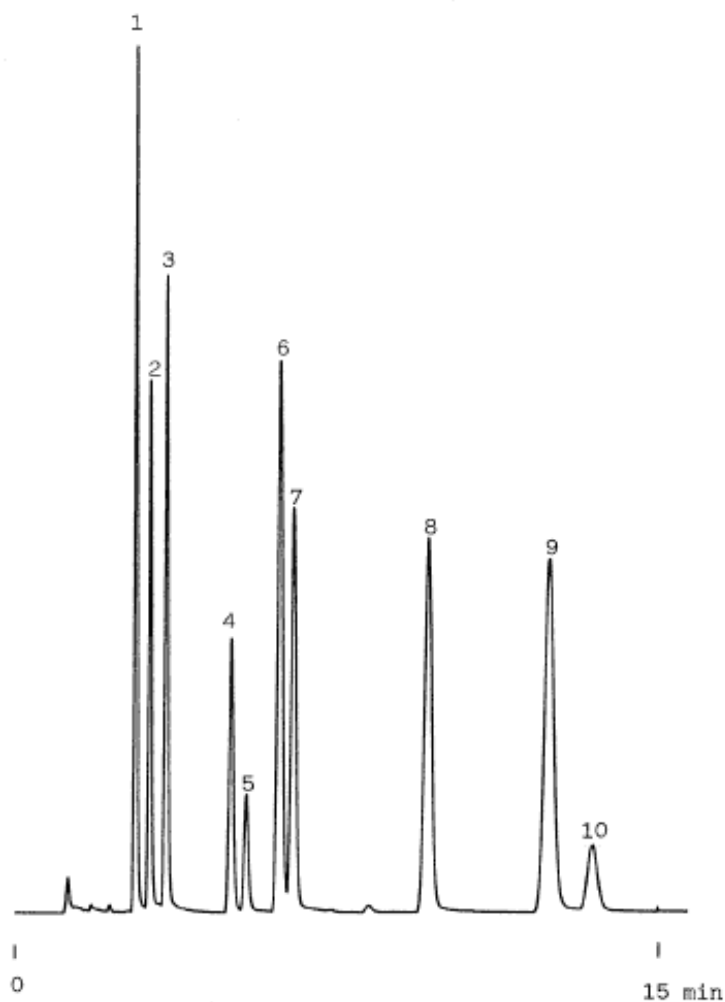

Column : Inertsil ODS-80A  
5 $\mu$ m  
150  $\times$  4.6mm I.D.

Column Cat. No. : 5020-01601

Eluent : CH<sub>3</sub>CN/0.1% H<sub>3</sub>PO<sub>4</sub>=60/40

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. MCPA  
2. Triclopyr  
3. MCPP  
4. Pyridaphenthion  
5. Methyl Dymron  
6. Mepronil  
7. Chloroneb  
8. Echlomezol  
9. Pencycuron  
10. Terbutcarb