


## Pesticides

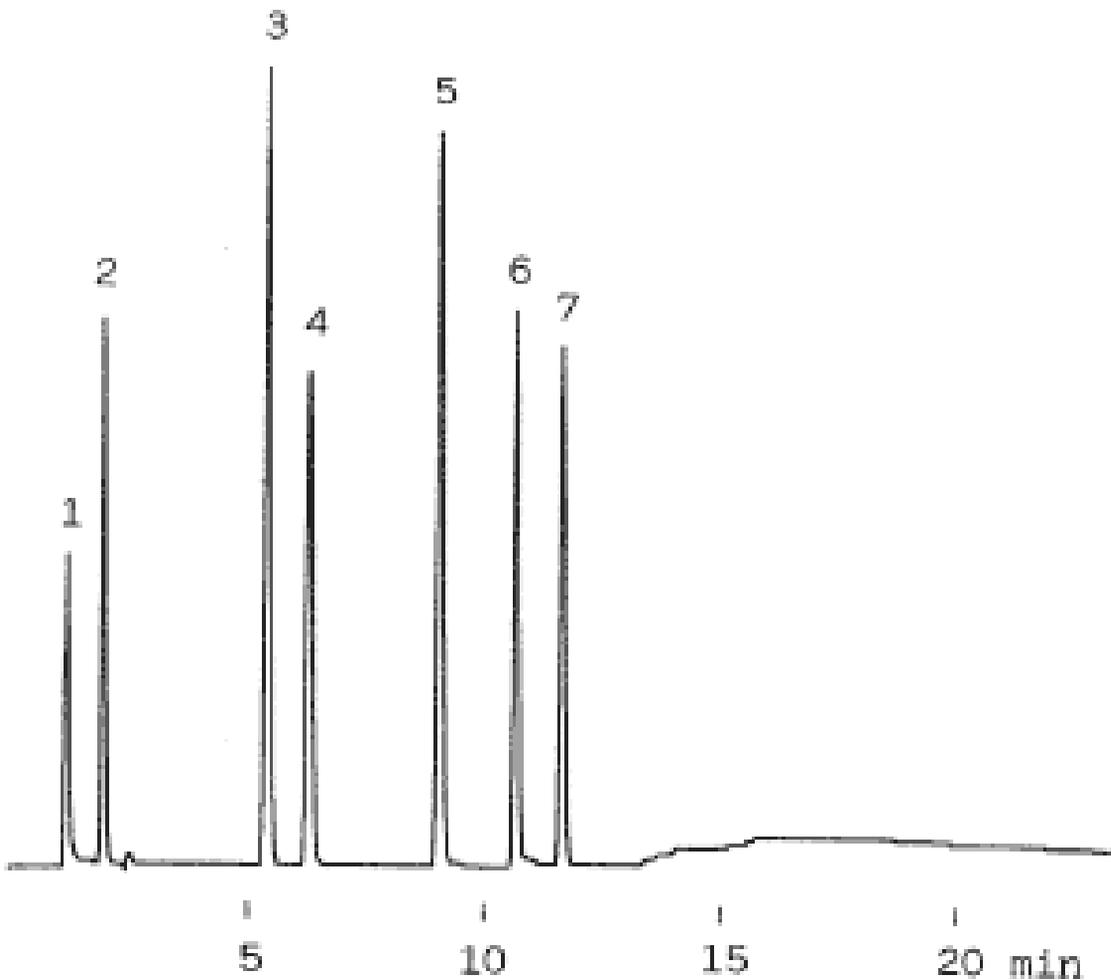

Column	: Inertsil ODS-3 5µm 150 × 4.6mmI.D.	Solutes	: 1. Asulam 2. Oxine-copper 3. Thiuram 4. MCPP 5. Iprodione 6. Bensulide 7. Pencycuron
Column Cat. No.	: 5020-13001		
Eluent	: A:0.2% H <sub>3</sub> PO <sub>4</sub> B:0.2% H <sub>3</sub> PO <sub>4</sub> in CH <sub>3</sub> CN A/B 50/50→(10min)→10/90 →(10min)→10/90		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 230 nm		
Injection Vol.	:		
Sample	:		