


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA290-0000

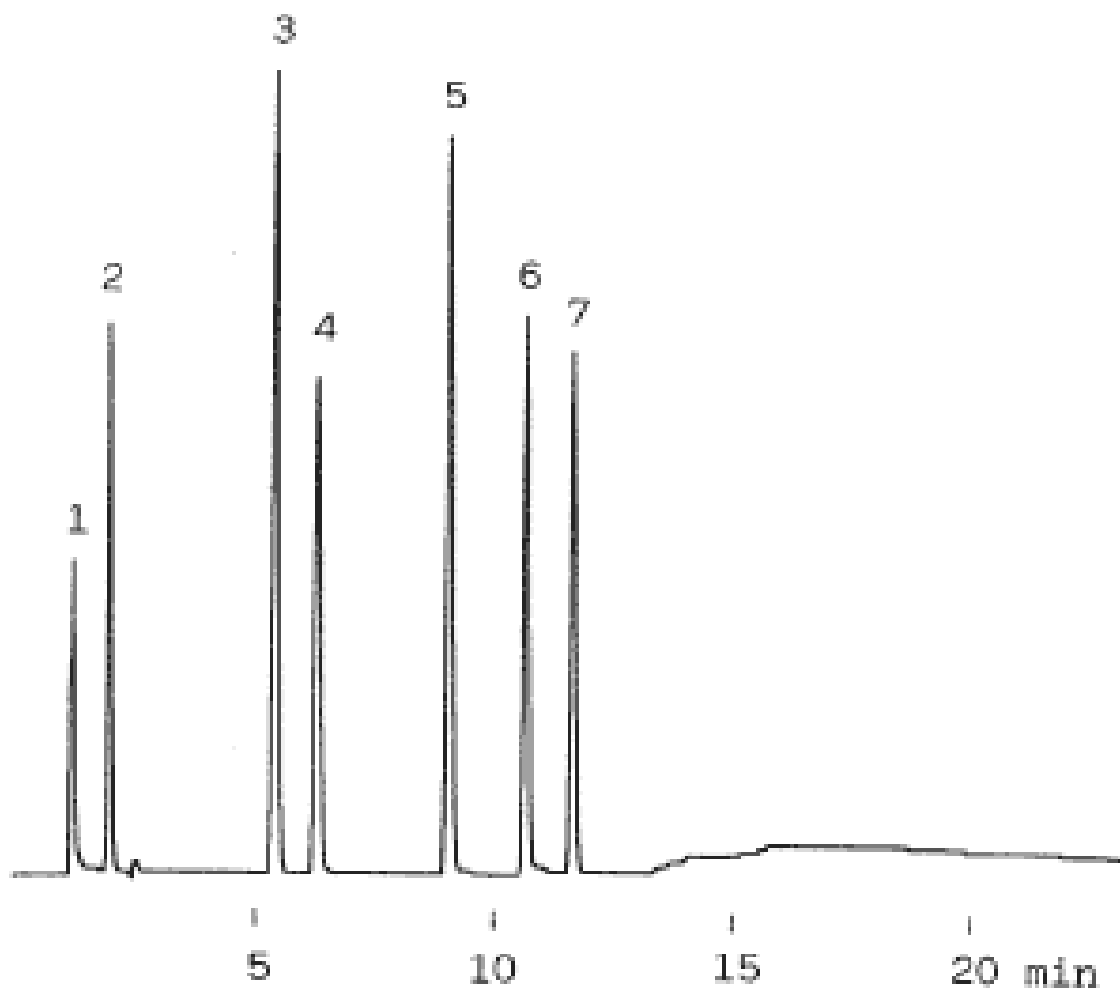

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-13001

Eluent : A: 0.2% H₃PO₄
B: 0.2% H₃PO₄ in CH₃CN
A/B 50/50 \rightarrow (10min) \rightarrow 10/90
 \rightarrow (10min) \rightarrow 10/90

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. Asulam
2. Oxine-copper
3. Thiuram
4. MCP
5. Iprodione
6. Bensulide
7. Pencycuron