


Analysis of Thiram

Data No. LA294-0000

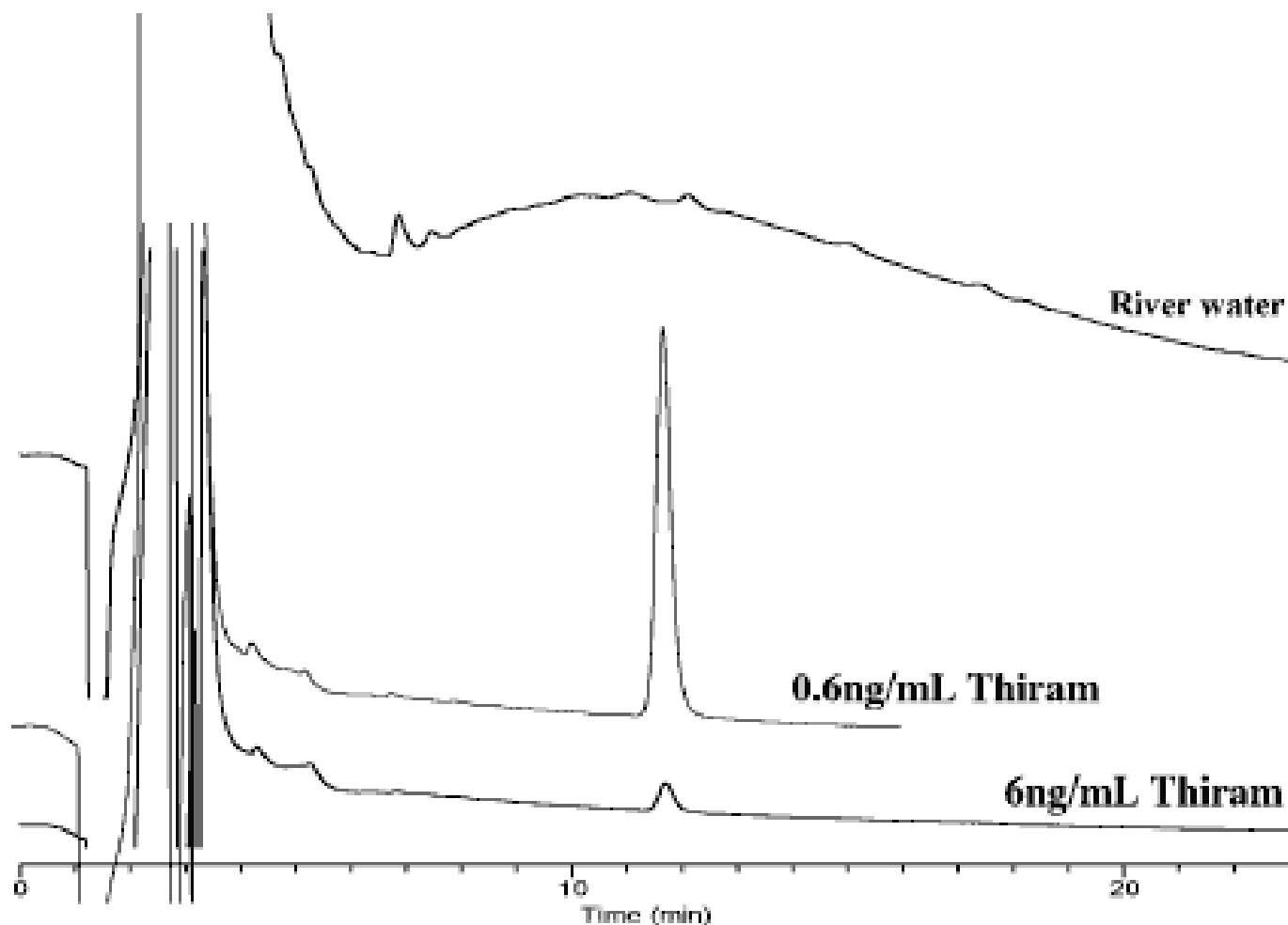

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D. Solutes : 1. Thiram (10 μ g/mL)

Column Cat. No. : 5020-01731

Eluent : Acetonitrile/50mM phosphate buffer (pH3.5) = 40/60

Flow rate : 1mL/min

Col. Temp. : 30°C

Detection : ECD 850mV vs. Ag/AgCl

Injection Vol. : 1 μ L

Sample : Pesticides