


Analysis of Thiram

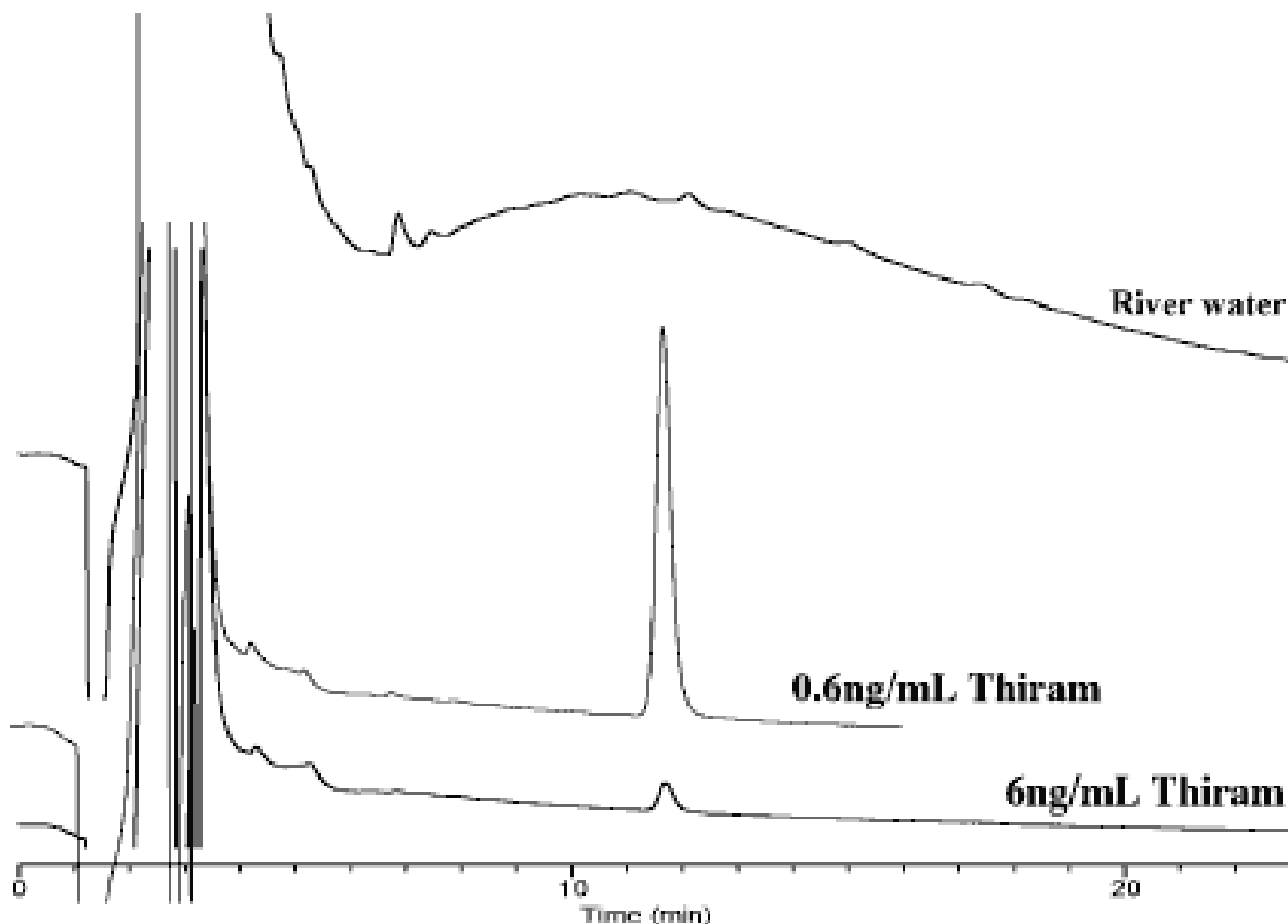

Column	: Inertsil ODS-3 5 μ m 150 \times 4.6mmI.D.	Solutes	: 1. Thiram (10 μ g/mL)
Column Cat. No. :	5020-01731		
Eluent	: Acetonitrile/50mM phosphate buffer (pH3.5) = 40/60		
Flow rate	: 1mL/min		
Col. Temp.	: 30°C		
Detection	: ECD 850mV vs. Ag/AgCl		
Injection Vol.	: 1 μ L		
Sample	: Pesticides		