


Analysis of Carbmates

Data No. LA296-0000

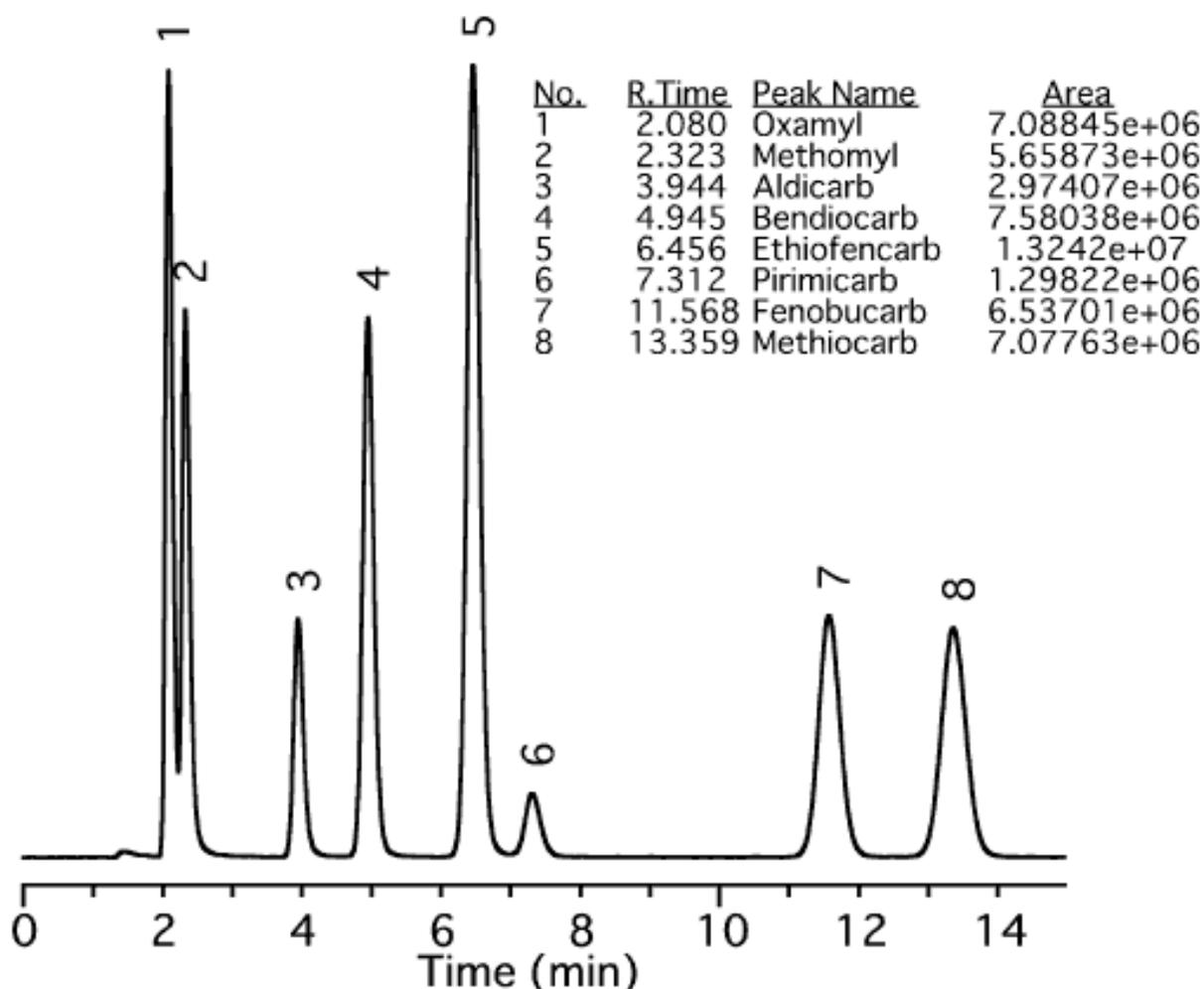

Column : Inertsil ODS-3
5µm
150 × 4.6mmI.D.

Column Cat. No. :

Eluent : 53% MeOH

Flow rate : 1mL/min

Col. Temp. : 40°C

Detection : FL Ex 345nm Em 450nm

POST REACTION NaOH 100? to OPA 50?

Injection Vol. : 25µL

Sample : Pesticides

Solutes : 1 . Oxamyl 2.83mg
2 . Methomyl 2.5mg
3 . Aldicarb 1.65mg
4 . Bendiocarb 2.5mg
5 . Ethiofencarb 5.10mg
6 . Pirimicarb 7.8mg
7 . Fenobucarb 2.0mg
8 . Methiocarb 2.4mg
in 50mL 52% MeOH