


Chlorobenzene

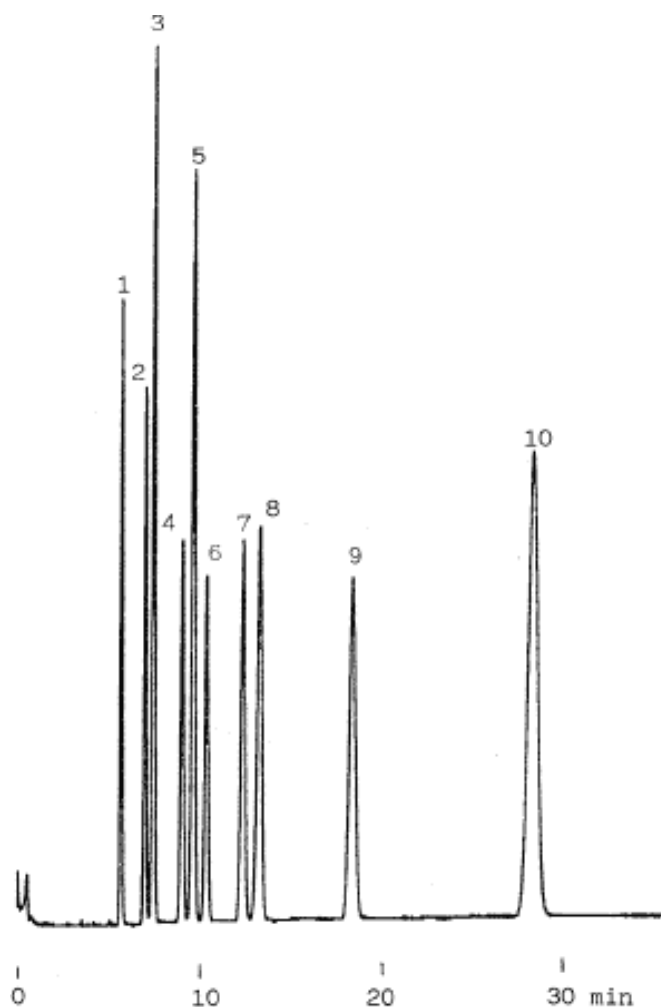

Column : Inertsil Ph
5 μ m
150 $\times$ 4.6mm I.D.

Column Cat. No. : 5020-01301

Eluent : CH₃CN/H₂O=50/50

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes	:	1. Chlorobenzene
		2. 1,2-Dichlorobenzene
		3. 1,4-Dichlorobenzene
		4. 1,3-Dichlorobenzene
		5. 1,2,3-Trichlorobenzene
		6. 1,3,5-Trichlorobenzene
		7. 1,2,3,4-Tetrachlorobenzene
		8. 1,2,4,5-Tetrachlorobenzene
		+1,2,3,5-Tetrachlorobenzene
		9. Pentachlorobenzene
		10. Hexachlorobenzene