
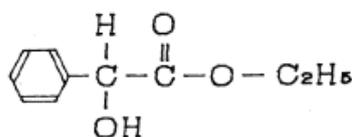


Mandelic acid ethyl ester

Data No. LA368-0000

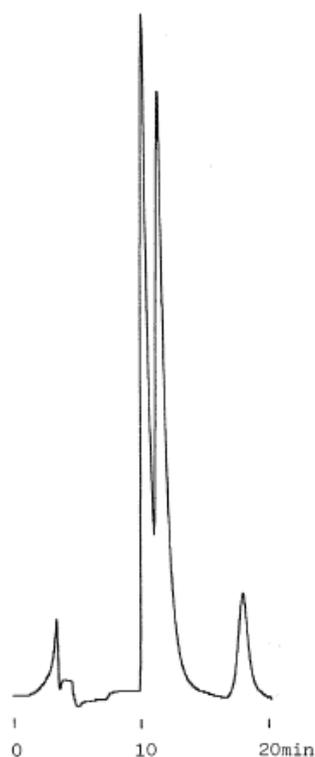

Column : Bioptic AV-1

5µm

150 × 4.6 mmI.D.

Column Cat. No. : 5020-02511

Eluent : 0.1M KH₂PO₄ +0.1% H₃PO₄

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 220nm

Injection Vol. : 1µL

Sample : Mandelic acid ethyl ester 10mg/3mL

Solutes : Mandelic acid ethyl ester