


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Fenoprofen

Data No. LA393-0000

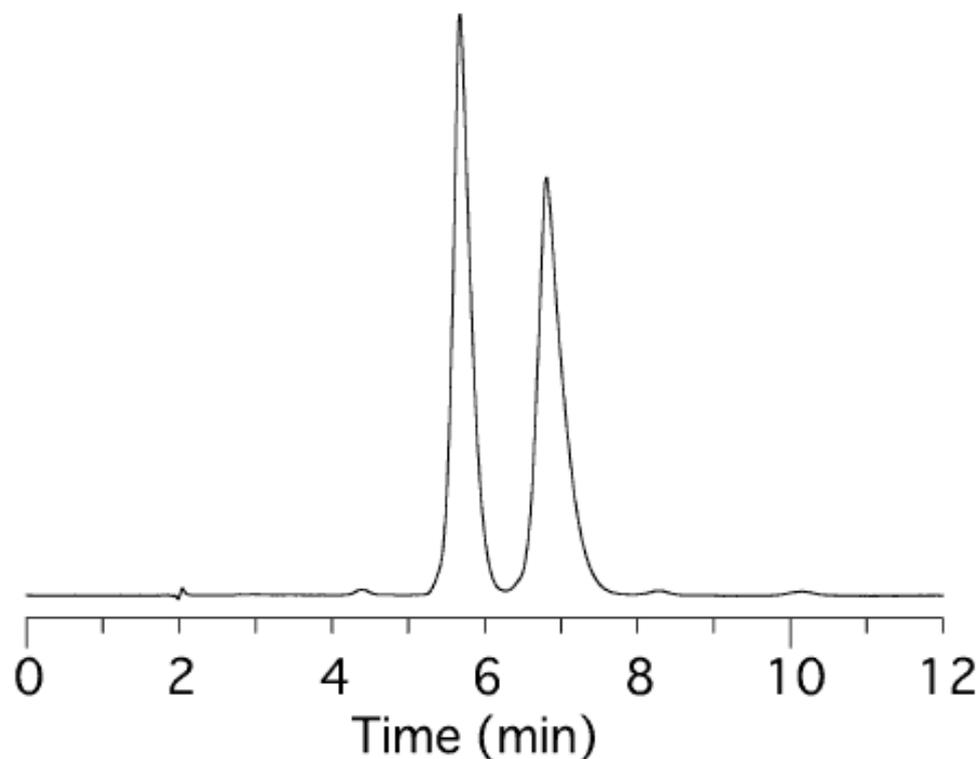

Column : Bioptic AV-1  
5 $\mu$ m  
150  $\times$  4.6 mm I.D. Solute : D,L-fenoprofen

Column Cat. No. : 5020-02511

Eluent : MeCN / 50mM phosphate buffer (pH7.3) = 5 / 95 (w / w)

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : D,L-fenoprofen