


DNBC-Alcohols

Data No. LA404-0000

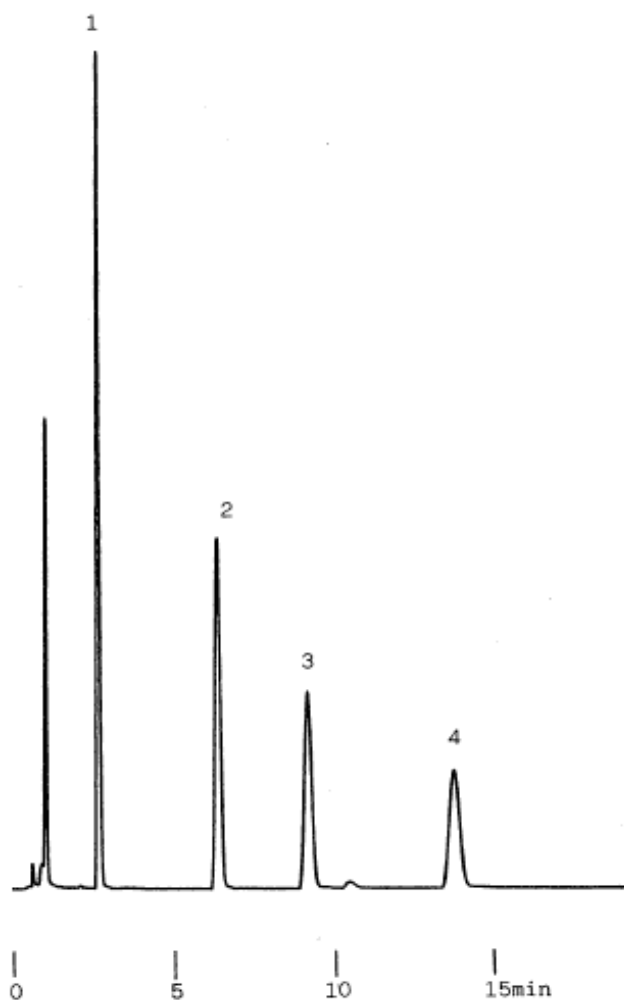

Column : Inertsil ODS-2
5 μ m
150 \times 4.6 mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/H₂O = 55/45

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : DNBC-Alcohol

Solutes : 1. DNBC-Ethylene glycol
2. DNBC-Ethanol
3. DNBC-n-Propanol
4. DNBC-n-Butanol